

KneeVision



Team S.T.A.R Labs

Agenda

01



SPRINT 2 RECAP

- Research Team
- Improvements
- Project Description
- Team Working Agreement
- Personas

02



PROJECT OVERVIEW

- MVP
- Languages and Tools
- Algorithms
- Diagrams

03



PRODUCT BACKLOG

- User Stories
- Sprint 3 Stories
- Test Cases

04



TEAM METRICS

- Team Velocity
- Burndown Charts
- Completed/Committed Ratio

05



PRODUCT

- Sprint 3 Retrospective
- API
- Live Application Demo

Our Research Team



Daniel Fox
**Team Leader/
Developer**



Manohar Killamsetti
**Scrum
Master/Developer**



Siming Li
Developer/Tester



Sarvesh Shah
Cloud Developer



Swatej Goud
Backend Developer



Deep Patel
ML Developer

Improvements

- Updated product backlog story points to correctly reflect Fibonacci numbers (slide 37) 
- Added averages along with the numbers in the metrics section (slides 51 – 54) 

Problem Statement

Knee X-ray imaging is commonly used to evaluate osteoarthritis, yet interpreting radiographic signs of disease severity can be subjective and inconsistent due to subtle visual differences and reliance on manual assessment. While advanced imaging techniques such as MRI provide greater detail, they are expensive and less accessible. This creates a need for scalable and consistent analysis methods that can assist in assessing knee osteoarthritis using widely available X-ray images.

Project Definition

- A web-based platform that uses deep learning to analyze knee X-ray images for knee osteoarthritis severity assessment.
- Intended for individuals, researchers, and clinicians seeking consistent and scalable analysis of knee X-ray data.
- Leverages publicly available datasets and established convolutional neural network (CNN) models.
- Designed as an AI-assisted decision-support and research tool, not a diagnostic system.
- Benefits include improved consistency, scalability, and reduced subjectivity in manual X-ray interpretation.

Project Description

Project Name:	KneeVision
Team:	S.T.A.R. Labs
Project Description:	KneeVision is a web-based platform that uses machine learning to analyze knee X-ray images, enabling faster and more accessible diagnostic insights. For researchers and medical professionals who upload knee X-ray images the KneeVision app is a web-based application that uses machine learning to analyze knee X-rays and assist in detecting potential conditions that provides faster, data-driven insights to support diagnosis and decision-making unlike traditional workflows that rely solely on manual review by specialists and can be time-consuming our application improves accessibility and efficiency by enabling quick image analysis and centralized result management
Benefit Outcomes:	• Faster turnaround time for preliminary image analysis • Improved accessibility to diagnostic insights • Centralized storage of user-specific medical images • Enhanced workflow for reviewing and managing X-ray data • Scalable foundation for integrating advanced AI/ML models
Github Link:	https://github.com/htmw/STARLABS

Team Agreement: Collaboration Standards

Meeting Pattern

Synchronous sessions twice weekly (Sun/Tue).

Mandated Agenda

documents for every session.

Continuous async sync via WhatsApp/Discord.

Decision Making

Majority-vote consensus based on researched evidence. Inclusive voicing of concerns during all discussion phases.

Technical Workflow

No direct pushes to Main. Isolated feature branches only.

Peer Approval. Every PR requires at least one secondary reviewer.

Jira Hygiene. Mandatory linking of work to Jira tracker IDs.

Blockers. Resolve locally & document attempts before escalation.

Team Agreement: Definition of “Done”

Category	Standard for Completion
Execution	Code runs without syntactical or runtime errors.
Compliance	Work meets all Industry Standards.
Verification	Work reviewed by multiple members & verified on multiple systems.
Documentation	All necessary docs tagged or updated along with code push.
Testing	All pre-defined test scenarios must be passed.

Persona - 1: The Practitioner

ORTHOPEDIC SURGEON | URBAN TEACHING HOSPITAL

Dr. Sarah Chen, 42

An expert in orthopedic medicine with 15 years of experience. Sarah operates in high-volume clinics where diagnostic speed is paramount to patient throughput and surgical triage.

KEY GOALS

- Rapid knee X-ray triage
- Reduce initial assessment time
- Validation for borderline cases
- Fatigue from patient volume
- Reader reporting variability
- Limited case analysis time



Persona - 2: The Research Coordinator

CLINICAL RESEARCH COORDINATOR

Alex Rodriguez, 29

A data-centric professional managing multi-site studies. Alex focuses on standardization, data integrity, and inter-rater reliability across large datasets.

STRATEGIC GOALS

Standardize severity grading	Efficient batch processing
Structured data exporting	Ensure rater reliability

OPERATIONAL PAINS

Manual grading	Lengthy training cycles
inconsistency	Audit trail requirements
Historical data re-analysis	



Persona - 3: The Learner



3RD YEAR MEDICAL STUDENT

Maya Patel, 24

A digital native in her clinical rotations. Maya uses the platform as a pedagogical tool to build pattern recognition skills and diagnostic confidence before her residency.

LEARNING GOALS

- Identify radiographic signs
- Immediate feedback practice
- Grade/Clinical correlation


ACCESS OBSTACLES

- Limited diverse case access
- Assessment uncertainty
- Limited attending time

What will our MVP be?

- Web-based interface
- Image preprocessing pipeline
- Deep learning model
- Severity assessment output
- Basic result visualization
- Mobile-responsive UI

Project Design (Figma)

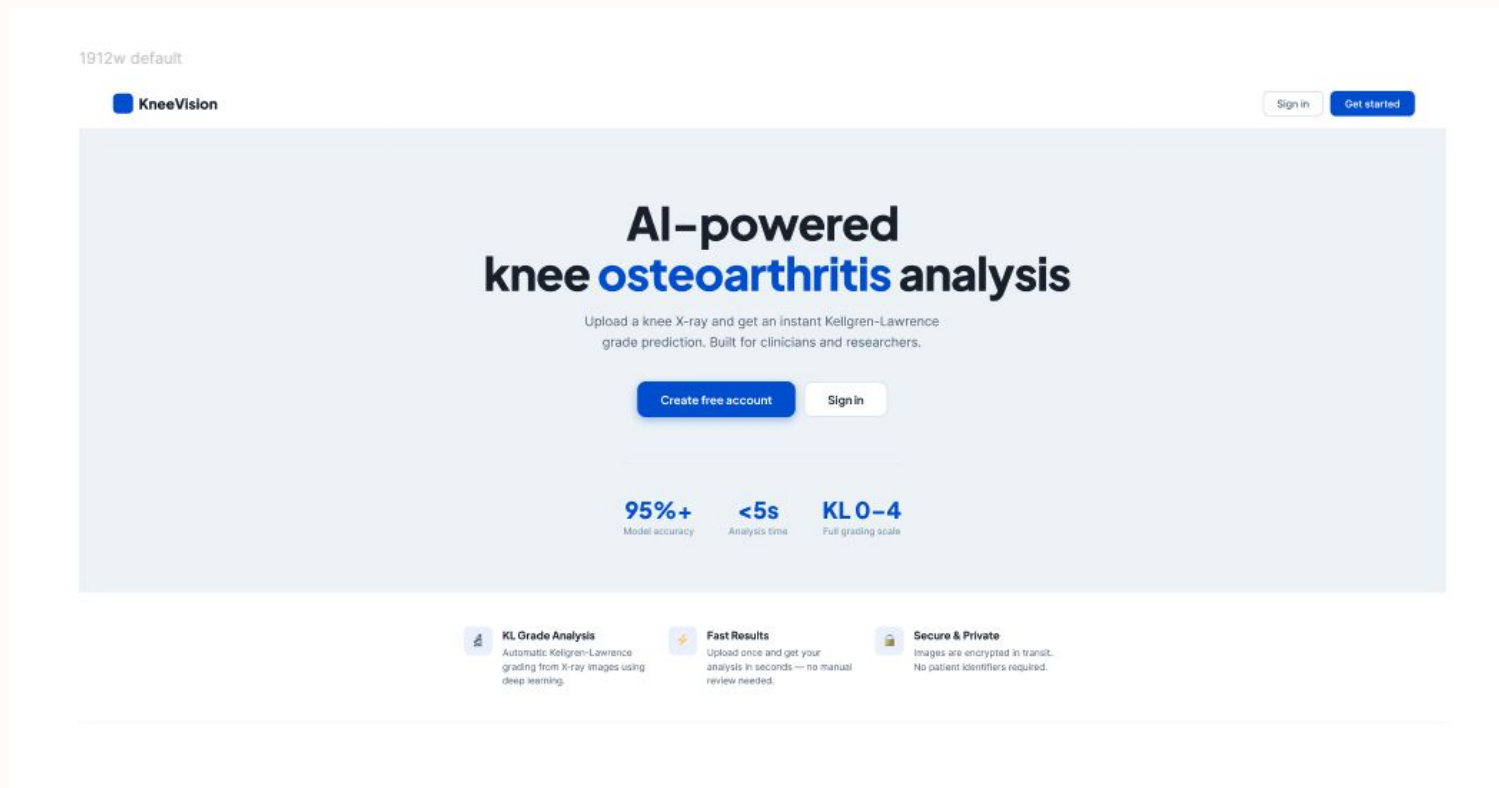

Tools & Technologies

Languages:  

Frontend:  

Backend:  

ML/DL Libraries: 

Database:  

Cloud: 

Tools:    

ML Inference: 

Misc. Libraries:  

FRONTEND

React · Vite · Tailwind CSS · TypeScript



React 18

Component-based UI

- Open-source JavaScript library for interactive UIs
- Reusable component architecture for maintainability
- State management with useState and useEffect hooks
- Pages: Landing, Login, Register, Upload, Gallery



Vite

Build toolchain

- Lightning-fast HMR and dev server startup
- Instant compilation during development
- Optimised production builds with tree-shaking
- Configured with Vite proxy to backend



Tailwind CSS

Utility-first styling

- Direct utility approach replaces custom CSS files
- Consistent design tokens across all components
- Responsive layouts with mobile-first breakpoints
- Dark-mode ready with zero extra config



TypeScript

Type-safe JavaScript

- Static typing catches errors at compile time
- Typed props and API response interfaces
- Improves IDE autocomplete and refactoring
- Compiled to JavaScript for browser execution

BACKEND

Node.js · Express.js · JWT · bcrypt



Node.js

Runtime environment

JavaScript runtime built on Chrome's V8 engine. Enables server-side JS execution with non-blocking I/O for high concurrency — perfect for handling multiple concurrent upload and analysis requests.



Express.js

Web framework

Minimal and flexible Node.js web framework. Handles all REST API routes: auth (register/login), image presign, raw file upload, image metadata registration, and image retrieval with JWT middleware.



JWT Auth

Stateless authentication

JSON Web Tokens secure every protected route via the requireAuth middleware. Tokens are signed with a secret, expire in 1 hour, and carry userId + email so any route knows who made the request.



bcrypt

Password hashing

Industry-standard adaptive hashing algorithm. Passwords are never stored in plain text — bcrypt hashes with salt rounds of 10, making brute-force attacks computationally infeasible.

DATABASE & CLOUD

MongoDB · AWS · Local filesystem



MongoDB

NoSQL document database

users collection

Stores email, hashed password, userId, timestamps

images collection

fileUrl, originalName, contentType, uploadedBy { userId, email }

predictions collection

imageId FK, klGrade, confidence, probabilities, modelVersion

Indexing

Unique index on users.email prevents duplicate accounts



AWS / Cloud

Cloud infrastructure (planned)

Amazon S3

Object storage for uploaded X-ray images at scale

EC2 / Lambda

Scalable compute for Express backend and ML inference

CloudFront

CDN for fast global delivery of static assets

Current dev

Local filesystem /uploads/ used during development

DEVOPS & TOOLS

Docker · GitHub · Jira · VS Code · Jest · Vitest



Docker

Containerises the Express backend and ML service so the environment is identical across dev, staging, and production machines.



GitHub

Version control and collaboration. Feature branches, pull requests, and code review. CI/CD pipeline triggers test suite on every push.



Jira

Agile project management. SCRUM board tracks sprints (SCRUM-18 register, SCRUM-19 login, SCRUM-20 image upload). Backlog prioritised by story points.



VS Code

Primary IDE with TypeScript, ESLint, and Prettier extensions. Integrated terminal runs both Vite dev server and Express backend simultaneously.



Jest + Vitest

Backend unit and integration tests with Jest (86 tests, ~89% coverage). Frontend component tests with Vitest + React Testing Library.



ESLint + Prettier

Code quality enforced at save. ESLint catches bugs, Prettier formats consistently. Husky pre-commit hooks block bad commits.

Model Pipeline Explanation

Dataset: Mendeley Dataset (8260 Images)

Model Plan:

Cleaning Dataset (Manually deleting the the defective images and cropping to ROI [Region of Interest])


Generating synthetic data using DCGAN (Deep Convolutional Generative Adversarial Network)

Preprocessing (Enhancing images using Adaptive Histogram Equalization [AHE])


Using Pre-Trained Models

Metrics Evaluation


Model Pipeline - Dataset

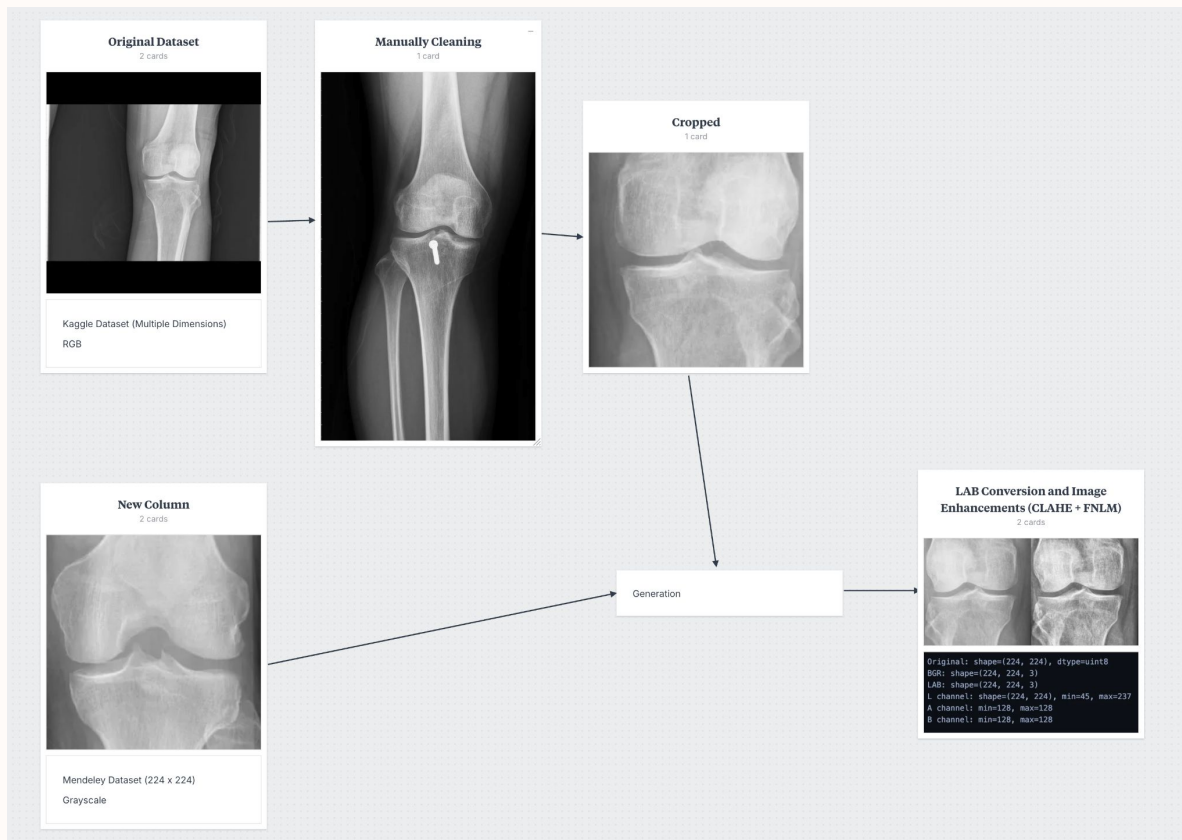

Generating Synthetic Data - DCGAN

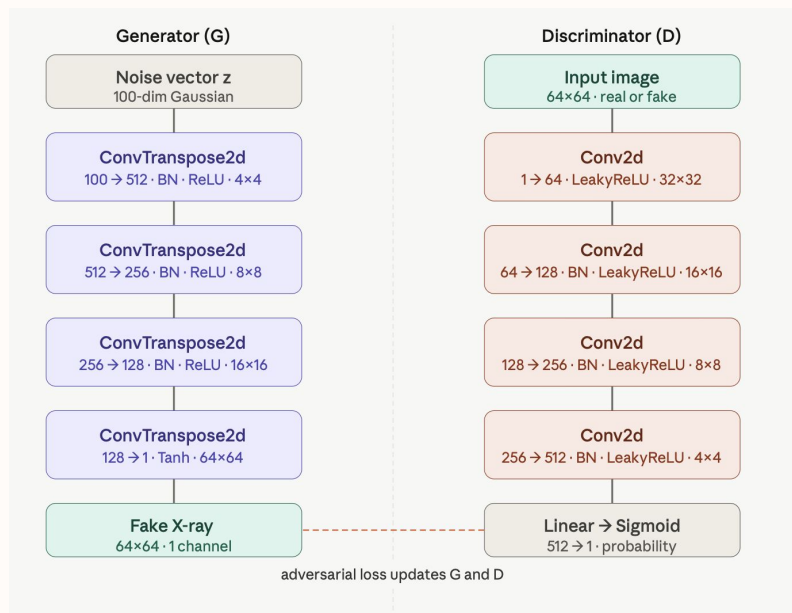


DCGAN Architecture

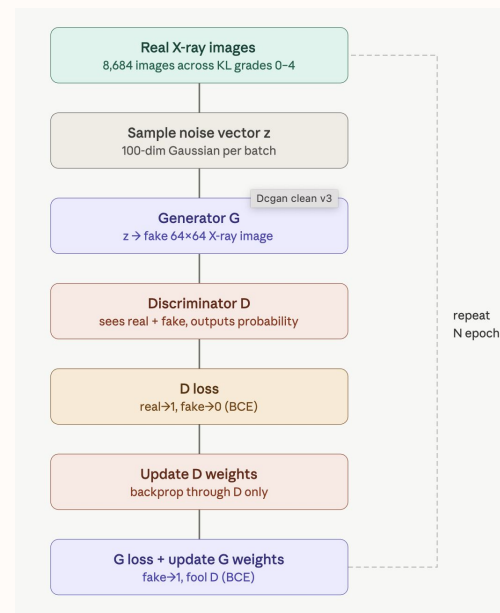


Generation Process

Generating Synthetic Data - DCGAN

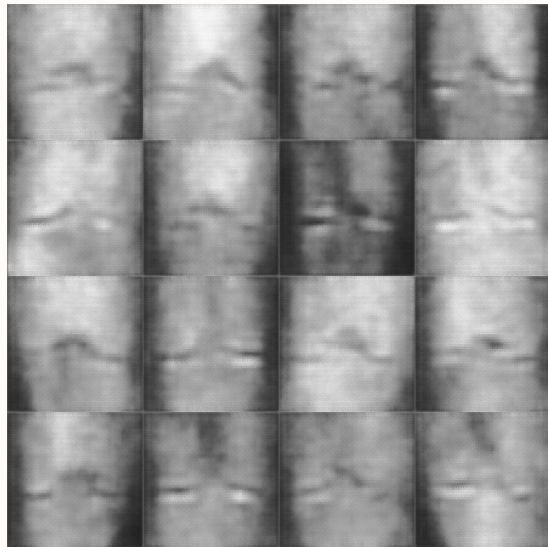

Epoch 0

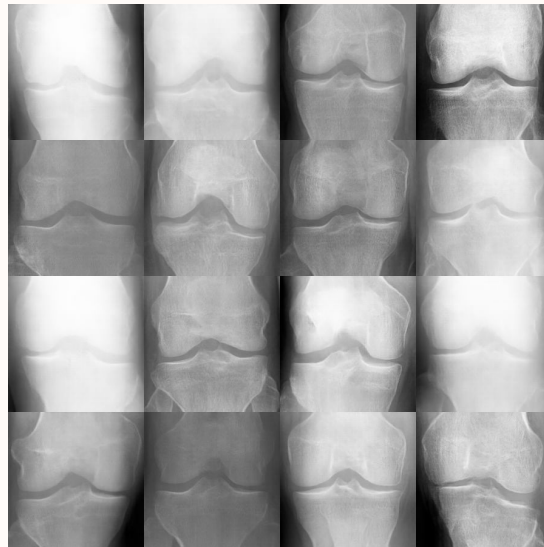



Epoch 200

Generating Synthetic Data - DCGAN




Model Overview

KneeVision — EfficientNet-B4 Training Process

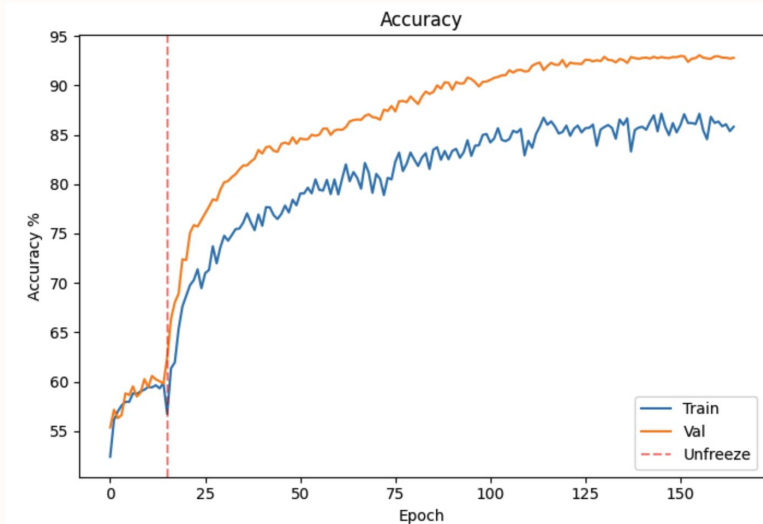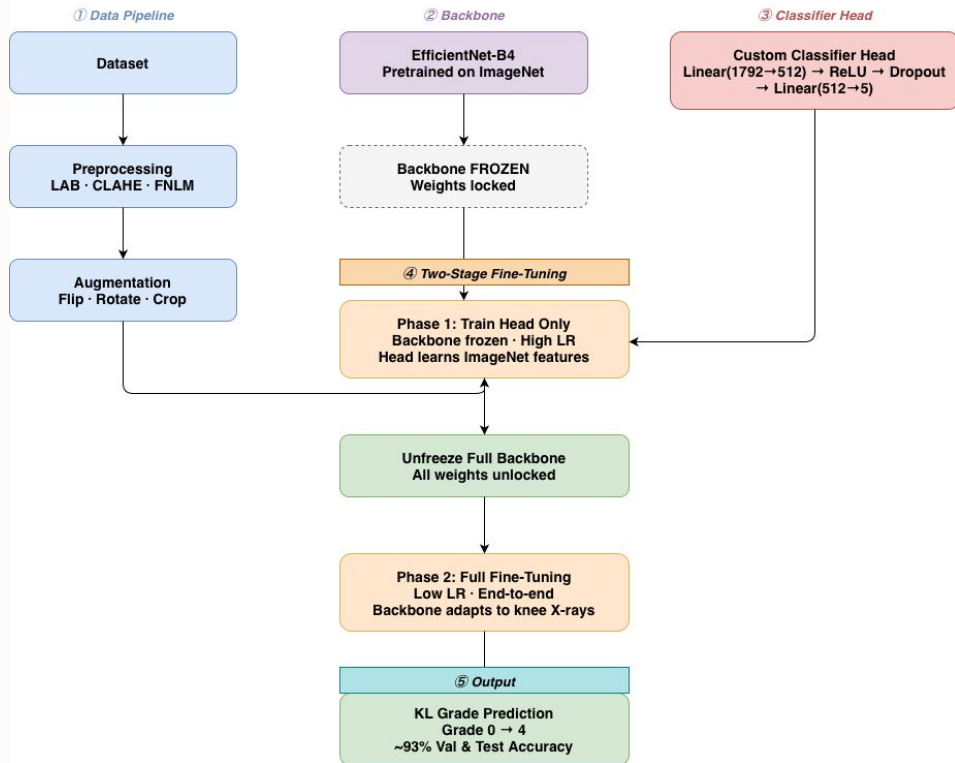


Context Diagram


Context Diagram

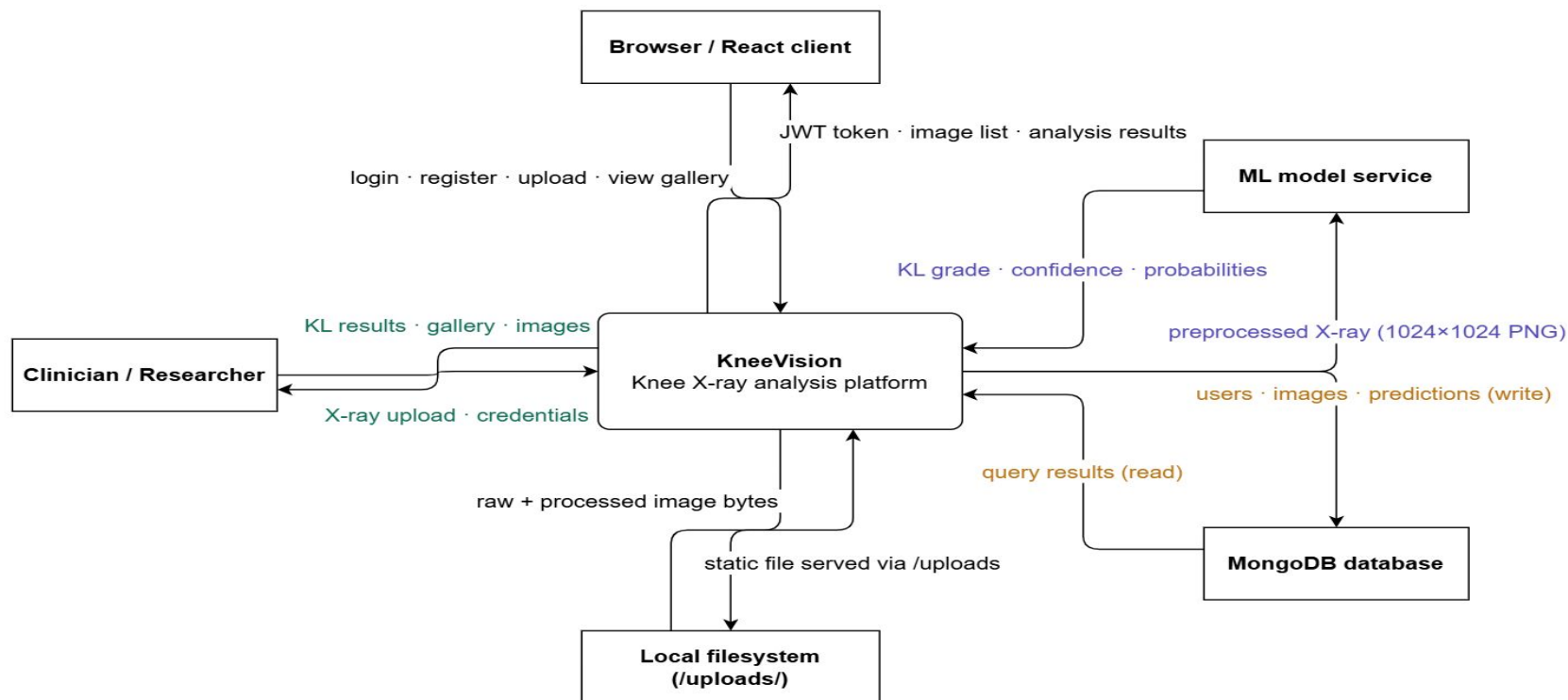

Class Diagram

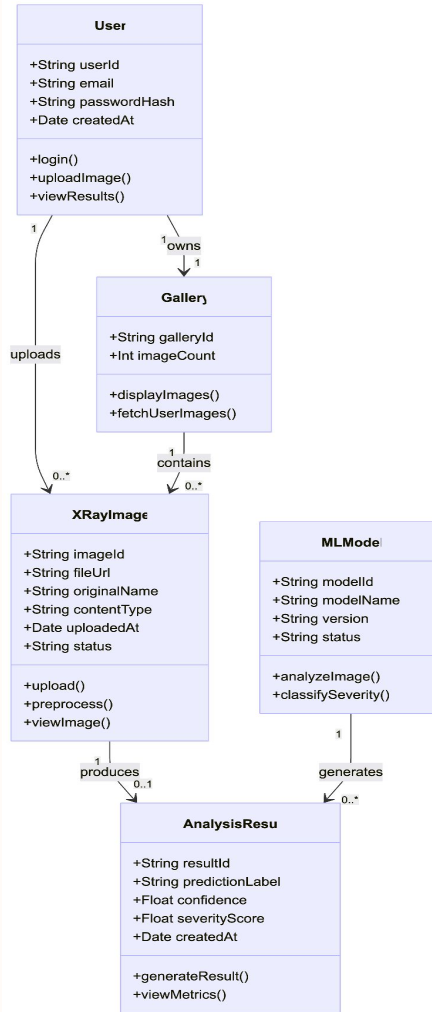


ER Diagram

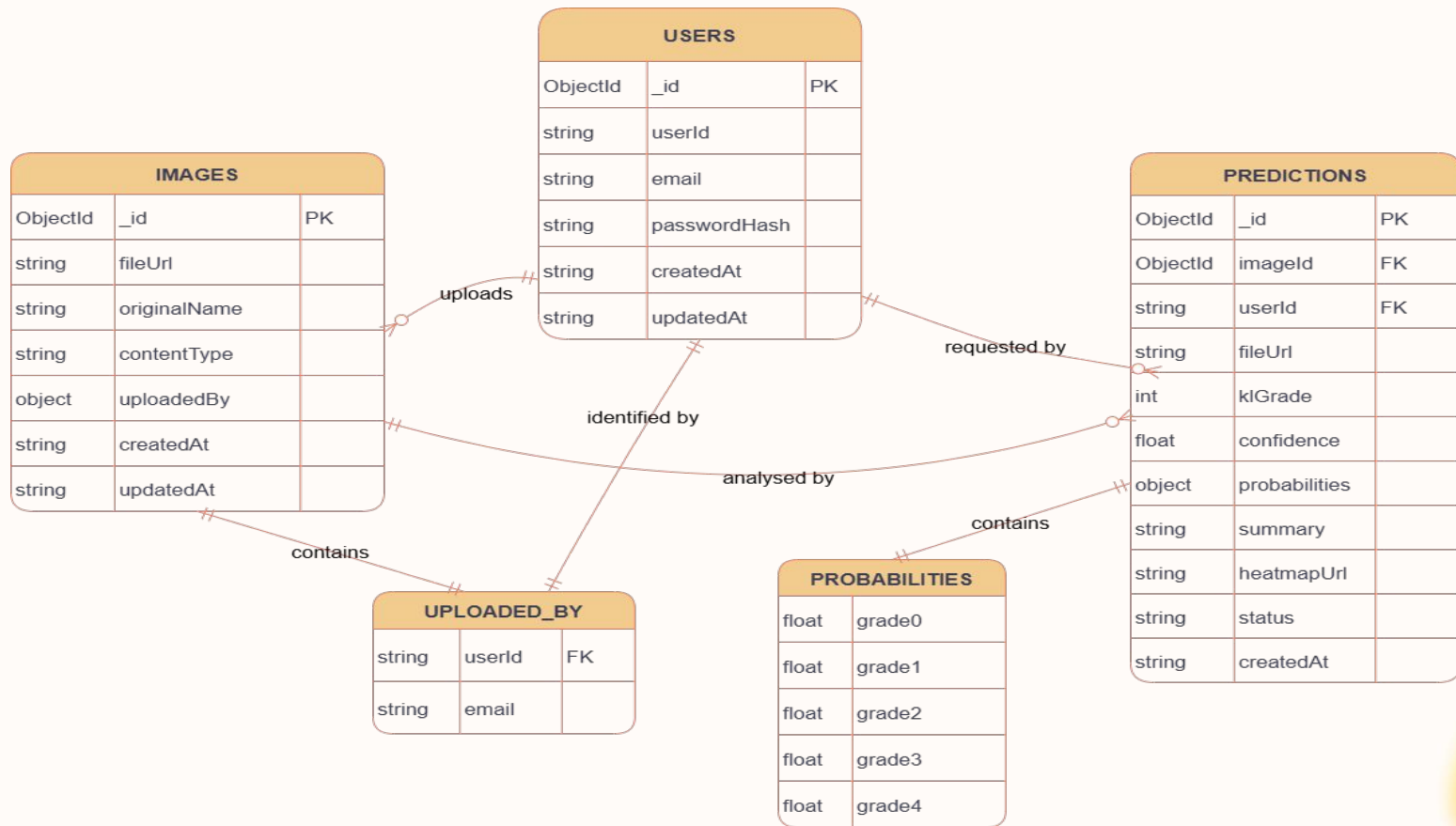

Architecture Diagram


KneeVision — System Design

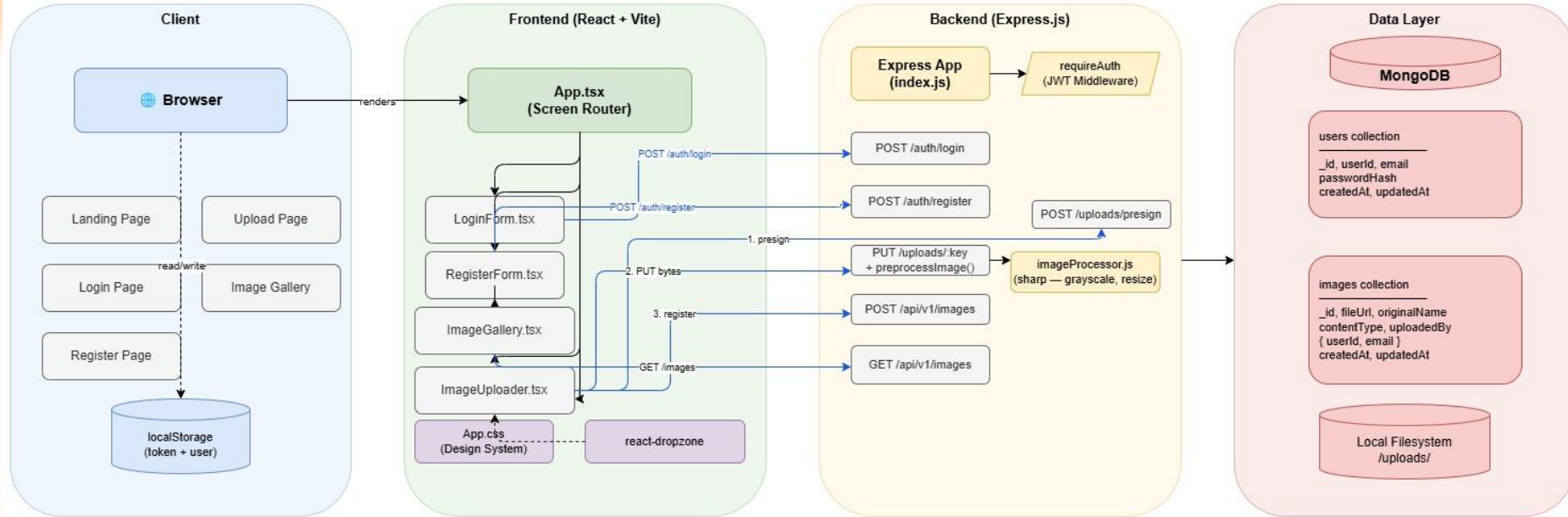


Legend

- HTTP request
- - - - - dependency

Sequence Diagram

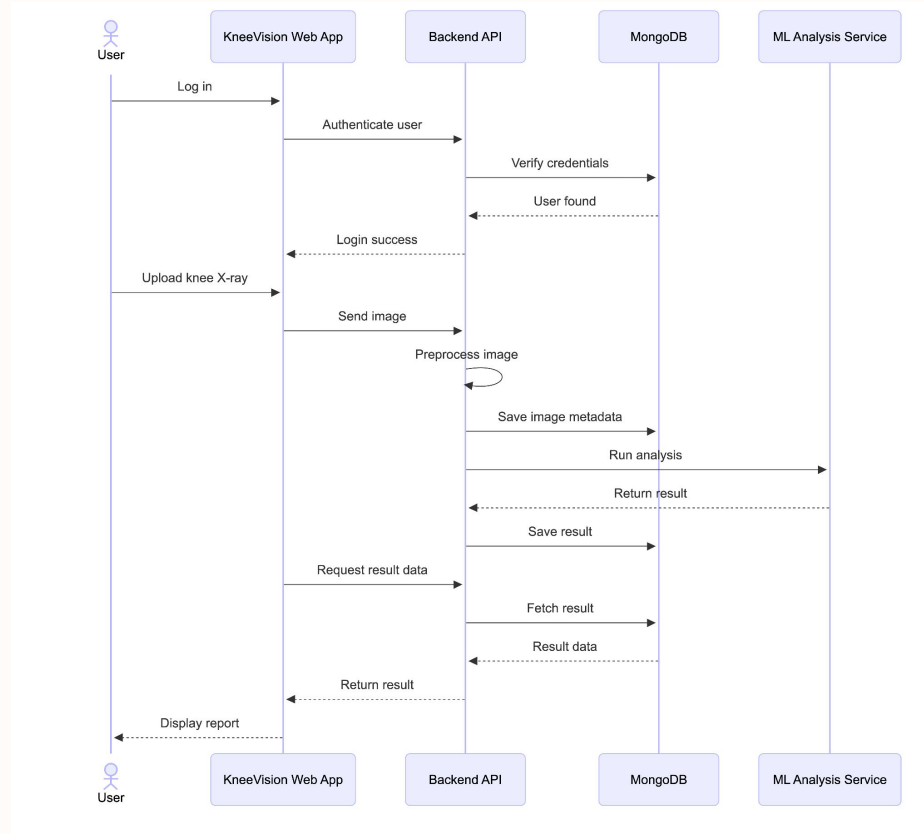


State Diagram

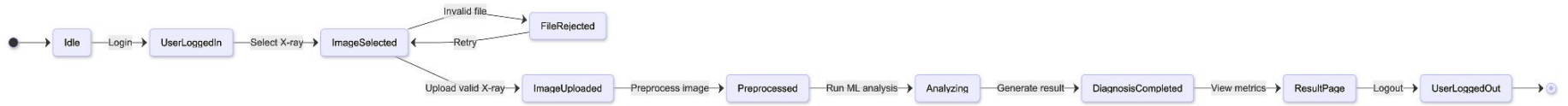


Sprint 2 Recap

- Integrated ML model and built prediction pipeline (upload → predict → results)
- Developed prediction API and model serving endpoint
- Built results page with KL grade, confidence, and visualizations
- Implemented analysis history page with filter and sort functionality
- Stored prediction results in database for persistence
- Conducted integration and cross-browser testing (~89% test coverage)

Sprint 1 Summary - Burndown Chart

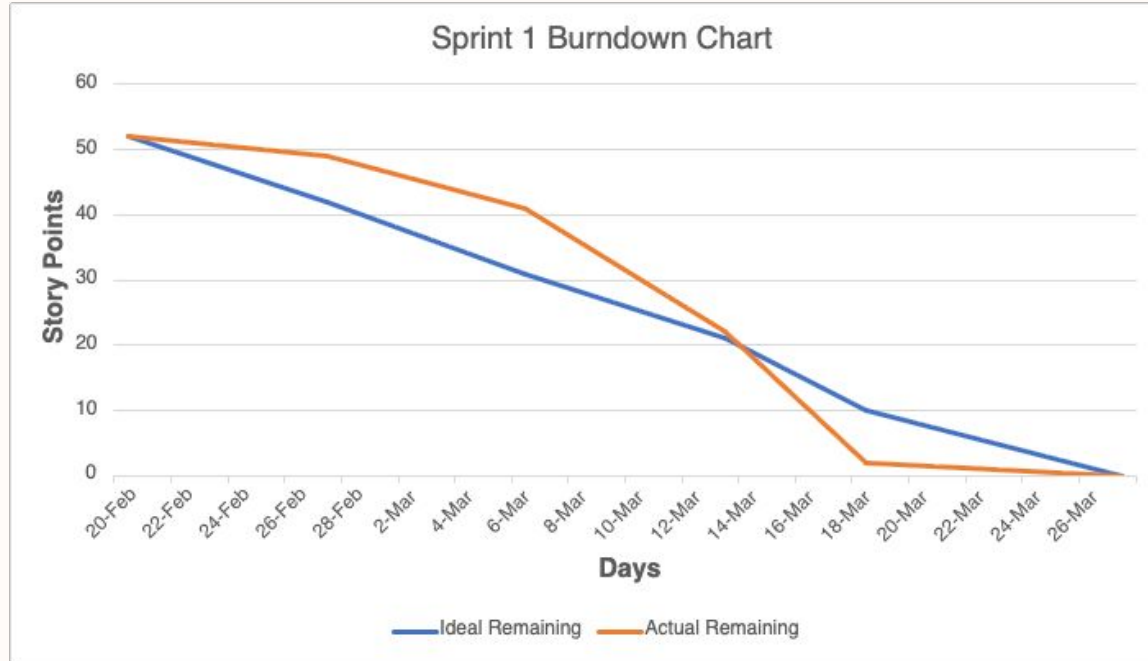

Sprint 2 Summary - Burndown Chart

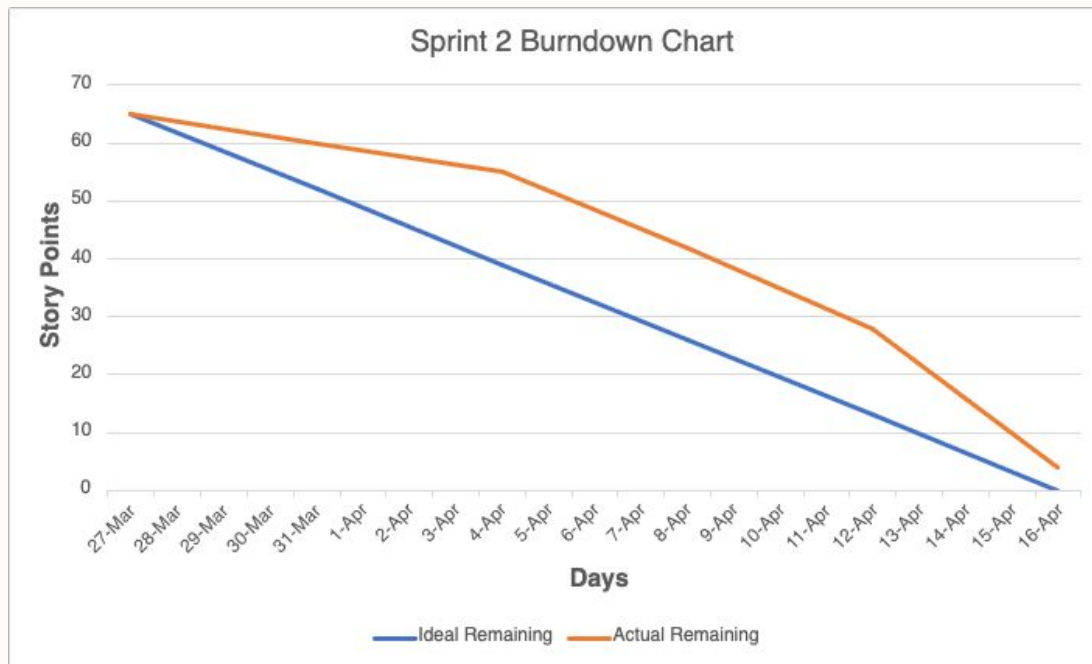

PRODUCT BACKLOG - USER STORIES

Sprint 1

Sprint 2

Sprint 3

US_ID	User Stories	Acceptance Criteria	Feature	Story Points
US_01	As an individual user, I want to upload knee X-ray images so that I can receive an AI-based assessment of osteoarthritis severity.	The patient can access an interface to upload an image. The patient can select and submit an image from device. The system accepts image files in specified formats: JPEG, JPG, PNG, DICOM. The uploaded image is stored and associated with the patient's case. A confirmation message is displayed upon successful upload.	Image Upload	13
US_02	As a research coordinator, I want to upload multiple X-ray images at once so that I can analyze datasets efficiently.	The patient can upload multiple images (up to 3) to one case. The system accepts multiple images in specified formats (JPEG, JPG, PNG). All uploaded images are associated with the same patient case. The patient receives a confirmation message after uploading each image. The system displays a preview of all uploaded images for review before analysis.	Image Upload	8

PRODUCT BACKLOG - USER STORIES

US_ID	User Stories	Acceptance Criteria	Feature	Story Points
US_03	As an individual user, I want to register an account so that I can store my basic information and access personalised services.	The user can access a registration page requiring: Full name, Email address, Password creation. Post-registration, the user can log in and update their profile information.	User Management	2
US_04	As an individual user, I want to log in to my account so that I can manage my cases and view diagnoses.	The user accesses a login page requiring: Email address, Password. Upon successful login, the user is directed to their dashboard.	User Management	1
US_05	As a user, I want to view my past cases and track their statuses so that I can monitor my health progress.	The user can access a 'Case History' section after logging in. All past cases are listed in chronological order (most recent first). Each case displays: case ID, case open date, status. The system updates case statuses in real-time. Clicking on a case entry navigates the user to that case's report.	User Management	2
US_06	As an individual user, I want to implement JWT authentication so that my session is secure.	JWT tokens are generated on login and validated on every protected route. Tokens expire after a defined period and require re-authentication. Invalid or expired tokens return a 401 Unauthorized response.	User Management	3

PRODUCT BACKLOG - USER STORIES

US_ID	User Stories	Acceptance Criteria	Feature	Story Points
US_07	As an individual user, I want to receive an OA severity diagnosis so that I know the knee condition I am facing.	After submitting the image, the user receives a diagnosis. The diagnosis is presented as the KL grade (0–4) and severity label. The system provides the diagnosis within a reasonable time (< 1 min). The diagnosis is stored and associated with the user's case.	AI Model Integration	21
US_08	As an individual user, I want to receive the most accurate possible diagnosis so that I can trust the assessment provided.	The system continuously improves the model and algorithms for better diagnostic results. The diagnosis accuracy is optimised based on model refinements and external validation. The confidence level of the diagnosis should be high (e.g., 95%).	AI Model Integration	13
US_09	As an individual user, I want to receive a brief description of the diagnosed condition so that I can understand the issue I am facing.	After receiving the diagnosis, the user is provided with a brief description of the diagnosed condition. The description is written in clear, understandable language. The description includes common symptoms and causes.	AI Model Integration	2
US_10	As an individual user, I want to receive treatment recommendations so that I can take action towards treating my knee condition.	After receiving the diagnosis, the user is provided with treatment recommendations. The description is written in clear, understandable language. The description includes common symptoms and causes of the condition.	AI Model Integration	3



Sprint 3 Backlog

US_ID	User Stories	Acceptance Criteria	Feature	Story Points
US_03	As an individual user, I want to register an account so that I can store my basic information and access personalised services.	The user can access a registration page requiring: Full name, Email address, Password creation. Post-registration, the user can log in and update their profile information.	User Management	2
US_04	As an individual user, I want to log in to my account so that I can manage my cases and view diagnoses.	The user accesses a login page requiring: Email address, Password. Upon successful login, the user is directed to their dashboard.	User Management	1
US_05	As a user, I want to view my past cases and track their statuses so that I can monitor my health progress.	The user can access a 'Case History' section after logging in. All past cases are listed in chronological order (most recent first). Each case displays: case ID, case open date, status. The system updates case statuses in real-time. Clicking on a case entry navigates the user to that case's report.	User Management	2
US_06	As an individual user, I want to implement JWT authentication so that my session is secure.	JWT tokens are generated on login and validated on every protected route. Tokens expire after a defined period and require re-authentication. Invalid or expired tokens return a 401 Unauthorized response.	User Management	3

Project Schedule - Sprint 3

Key	Issue Type	Summary	Status	Assignee	Sprint
SCRUM-59	Task	Add educational tooltips throughout app	Done	Swatej	Sprint 3
SCRUM-60	Task	Make UI responsive (mobile/tablet/desktop)	Done	Daniel	Sprint 3
SCRUM-33	Task	Deploy to development environment	Done	Daniel	Sprint 3
SCRUM-43	Task	Test prediction API (< 7 sec response)	Done	Siming	Sprint 3
SCRUM-47	Task	Show confidence score with progress bar	Done	Deep	Sprint 3
SCRUM-61	Task	Add error logging (backend + frontend)	Done	Swatej	Sprint 3
SCRUM-62	Task	Conduct UAT with 3-5 external testers	Done	Manohar	Sprint 3
SCRUM-64	Task	Performance testing (50 concurrent users)	Done	Manohar	Sprint 3
SCRUM-65	Task	Security audit (OWASP top 10)	Done	Manohar	Sprint 3
SCRUM-66	Task	Fix all critical and high-priority bugs	Done	Daniel	Sprint 3
SCRUM-67	Task	Deploy frontend to Vercel/Netlify	Done	Manohar	Sprint 3
SCRUM-68	Task	Deploy backend to Render/Railway	Done	Manohar	Sprint 3
SCRUM-69	Task	Deploy ML model endpoint	Done	Manohar	Sprint 3

Project Schedule - Sprint 3

Key	Issue Type	Summary	Status	Assignee	Sprint
SCRUM-70	Task	Set up monitoring and error tracking	Done	Siming	Sprint 3
SCRUM-71	Task	Final end-to-end testing in production	Done	Daniel	Sprint 3
SCRUM-72	Task	Create privacy policy page	Done	Deep	Sprint 3
SCRUM-73	Task	Finalize technical paper (incorporate feedback)	Done	Deep	Sprint 3
SCRUM-74	Task	Write deployment manual (step-by-step)	Done	Manohar	Sprint 3
SCRUM-75	Task	Write user manual with screenshots	Done	Manohar	Sprint 3
SCRUM-76	Task	Record final MVP demo video (5-10 min)	Done	Sarvesh	Sprint 3
SCRUM-77	Task	Prepare final presentation slides	Done	Sarvesh	Sprint 3

Test Cases

Test Summary

Module	Total	Passed	Failed	Coverage
Backend — imageProcessor.js	11	11	0	~95%
Backend — db.js	6	6	0	~90%
Backend — API routes (index.js)	52	52	0	~90%
Frontend — LoginForm	10	10	0	~92%
Frontend — RegisterForm	10	10	0	~91%
Frontend — ImageGallery	8	8	0	~90%
Frontend — ImageUploader	8	8	0	~85%
Frontend — Pages	16	16	0	~88%
TOTAL	121	121	0	~89%

Stack==> Node.js | Express | MongoDB | React | TypeScript | Jest | Vitest

Test Cases - Backend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
1	Returns a Buffer	1. Create a 100×100 RGB PNG buffer using sharp. 2. Call preprocessImage(buffer). 3. Inspect the returned value.	preprocessImage() output is a Buffer instance	✓ PASS
2	Outputs PNG by default	1. Call preprocessImage(buffer) with no format option. 2. Inspect output metadata using sharp(result).metadata().	Default format produces image/png metadata	✓ PASS
3	Converts image to grayscale	1. Create a colour PNG buffer. 2. Call preprocessImage(buffer). 3. Read output metadata with sharp.	Output has 1 channel (grayscale)	✓ PASS
4	No upscale for small images	1. Call preprocessImage(buffer, { width: 1024, height: 1024 }). 2. Read output metadata.	50x50 image stays ≤50px when target is 1024x1024	✓ PASS
5	Resizes large image to fit target	1. Call preprocessImage(buffer, { width: 512, height: 512 }). 2. Read output metadata.	2000x2000 image resized to ≤512x512 when target is 512	✓ PASS
6	Preserves aspect ratio	1. Call preprocessImage(buffer, { width: 200, height: 200 }). 2. Compute width/height ratio of output.	200x100 (2:1) input produces ~2:1 ratio output	✓ PASS
7	Outputs JPEG when format=jpeg	1. Call preprocessImage(buffer, { format: 'jpeg' }). 2. Read output metadata.	format:'jpeg' option produces jpeg metadata	✓ PASS
8	Outputs WebP when format=webp	1. Call preprocessImage(buffer, { format: 'webp' }). 2. Read output metadata.	format:'webp' option produces webp metadata	✓ PASS
9	Normalize option does not throw	1. Call preprocessImage(buffer, { normalize: true }). 2. Await the result.	normalize:true runs without error	✓ PASS
10	Throws on invalid buffer	1. Create Buffer.from('not-an-image'). 2. Call preprocessImage(badBuffer). 3. Await and expect rejection.	Non-image buffer rejects with 'Image preprocessing failed'	✓ PASS
11	Throws on null input	1. Call preprocessImage(null). 2. Await and expect rejection.	null input causes rejection	✓ PASS
12	Throws when MONGODB_URI missing	1. Delete process.env.MONGODB_URI. 2. Call connectToMongo(). 3. Await the promise.	connectToMongo() rejects with 'MONGODB_URI is not set'	✓ PASS
13	Connects and returns db instance	1. Set process.env.MONGODB_URI to a valid URI. 2. Call connectToMongo(). 3. Inspect return value.	Valid URI resolves to a db object	✓ PASS
14	Creates unique index on users.email	1. Call connectToMongo(). 2. Check that db.collection('users').createIndex was called.	createIndex called with {email:1}, {unique:true}	✓ PASS
15	Caches db on repeated calls	1. Call connectToMongo() twice. 2. Check client.connect call count.	client.connect called exactly once after two connectToMongo() calls	✓ PASS

Test Cases - Backend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
16	getDb throws before connect	1. Import getDb from fresh module. 2. Call getDb() immediately.	getDb() throws 'MongoDB not connected' before connectToMongo()	✓ PASS
17	getDb returns db after connect	1. Call connectToMongo(). 2. Call getDb(). 3. Inspect return value.	getDb() returns db object after successful connect	✓ PASS
18	GET /health returns 200	1. Send GET request to /health.	Response body is {status:'ok'}	✓ PASS
19	Register — valid input returns 201	1. POST /api/v1/auth/register with { email: 'test@example.com', password: 'password123' }. 2. Inspect response.	Returns id, email, createdAt on valid body	✓ PASS
20	Register — missing email → 400	1. POST /api/v1/auth/register with { password: 'password123' } only.	Returns 400 with 'required' message	✓ PASS
21	Register — missing password → 400	1. POST /api/v1/auth/register with { email: 'test@example.com' } only.	Returns 400	✓ PASS
22	Register — invalid email → 400	1. POST /api/v1/auth/register with { email: 'not-an-email', password: 'password123' }.	Returns 400 with 'Invalid email format'	✓ PASS
23	Register — short password → 400	1. POST /api/v1/auth/register with { email: 'test@example.com', password: 'short' }.	Returns 400 with '8 characters' message	✓ PASS
24	Register — duplicate user → 409	1. POST /api/v1/auth/register with existing email.	Returns 409 with 'User already exists'	✓ PASS
25	Register — email normalized	1. POST with { email: 'USER@EXAMPLE.COM', password: 'password123' }. 2. Check response body email field.	USER@EXAMPLE.COM stored as user@example.com	✓ PASS
26	Login — valid credentials → 200	1. POST /api/v1/auth/login with correct email and password.	Returns JWT token and user object	✓ PASS
27	Login — empty body → 400	1. POST /api/v1/auth/login with empty body {}.	Returns 400	✓ PASS
28	Login — unknown user → 401	1. POST /api/v1/auth/login with unregistered email.	Returns 401 with 'Invalid email or password'	✓ PASS
29	Login — wrong password → 401	1. POST /api/v1/auth/login with correct email but wrong password.	Returns 401	✓ PASS
30	Presign — valid input → 200	1. POST /api/v1/uploads/presign with { filename: 'xray.png', contentType: 'image/png' }.	Returns uploadUrl, fileUrl, method:PUT	✓ PASS

Test Cases - Backend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
31	Presign — missing filename → 400	1. POST /api/v1/uploads/presign with { contentType: 'image/png' } only.	Returns 400	✓ PASS
32	Presign — missing contentType → 400	1. POST /api/v1/uploads/presign with { filename: 'xray.png' } only.	Returns 400	✓ PASS
33	Images — valid input → 201	1. POST /api/v1/images with { fileUrl: '/uploads/test.png', originalName: 'xray.png', contentType: 'image/png' }.	Returns id, fileUrl, originalName, contentType	✓ PASS
34	Images — missing fileUrl → 400	1. POST /api/v1/images with { originalName: 'xray.png' } only.	Returns 400 with 'fileUrl is required'	✓ PASS
35	GET /images — no token → 401	1. Send GET /api/v1/images with no Authorization header.	Returns 401 with 'Authorization token is required'	✓ PASS
36	GET /images — malformed token → 401	1. Send GET /api/v1/images with Authorization: Bearer not.a.real.token.	Returns 401 with 'Invalid or expired token'	✓ PASS
37	GET /images — empty gallery → 200 []	1. mockFind returns []. 2. Send GET /api/v1/images with valid Bearer token.	Returns 200 with empty array []	✓ PASS
38	GET /images — returns images mapped correctly	1. mockFind returns 2 fake image docs. 2. Send GET /api/v1/images with valid Bearer token.	Returns 200, body length 2, first item has id 'img001' and fileUrl '/uploads/scan.png'	✓ PASS
39	GET /images — DB error → 500	1. mockFind rejects with Error('DB error'). 2. Send GET /api/v1/images with valid Bearer token.	Returns 500 with 'Failed to fetch images'	✓ PASS
40	POST /predict — no token → 401	1. POST /api/v1/predict with body { fileUrl: '/uploads/scan.png' } and no token.	Returns 401	✓ PASS
41	POST /predict — missing fileUrl → 400	1. POST /api/v1/predict with valid token and empty body {}.	Returns 400 with 'fileUrl is required'	✓ PASS
42	POST /predict — invalid fileUrl prefix → 400	1. POST /api/v1/predict with { fileUrl: '/etc/passwd' }.	Returns 400 with 'Invalid fileUrl'	✓ PASS
43	POST /predict — file not on disk → 404	1. mockFsAccess rejects (ENOENT). 2. POST /api/v1/predict with { fileUrl: '/uploads/scan.png' }.	Returns 404 with 'Uploaded file not found'	✓ PASS
44	POST /predict — image not in DB → 404	1. mockFsAccess resolves. mockCallMIPredict resolves. mockFindOne returns null. 2. POST /api/v1/predict.	Returns 404 with 'Image not found'	✓ PASS
45	POST /predict — returns cached prediction	1. mockFindOne returns imageDoc then existing prediction. 2. POST /api/v1/predict. 3. Check mockInsertOne was not called.	Returns 200, id='pred999', klGrade=1. mockInsertOne called 0 times	✓ PASS

Test Cases - Backend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
46	POST /predict — saves new prediction and returns ML result	1. mockFindOne returns imageDoc then null. 2. mockCallMIPredict returns {klGrade:2, confidence:0.87, label:'Moderate'}. 3. mockInsertOne returns insertedId 'pred001'. 4. POST /api/v1/predict.	Returns 200 with {id:'pred001', imageId:'img001', klGrade:2, confidence:0.87}. mockInsertOne called once	✓ PASS
47	POST /predict — ML service error → 500	1. mockFsAccess resolves. 2. mockCallMIPredict throws Error('ML service down'). 3. POST /api/v1/predict.	Returns 500 with 'Prediction failed'	✓ PASS
48	GET /predictions — no token → 401	1. Send GET /api/v1/predictions with no Authorization header.	Returns 401	✓ PASS
49	GET /predictions — invalid token → 401	1. Send GET /api/v1/predictions with Authorization: Bearer invalid.token.here.	Returns 401	✓ PASS
50	GET /predictions — empty list → 200 []	1. mockFind returns []. 2. Send GET /api/v1/predictions with valid Bearer token.	Returns 200 with empty array []	✓ PASS
51	GET /predictions — returns predictions mapped correctly	1. mockFind returns 2 fake prediction docs. 2. Send GET /api/v1/predictions with valid Bearer token.	Returns 200, body length 2. First item: id='pred001', result.klGrade=2. Second item: result.klGrade=0	✓ PASS
52	GET /predictions — DB error → 500	1. mockFind rejects with Error('DB connection lost'). 2. Send GET /api/v1/predictions with valid Bearer token.	Returns 500 with 'Failed to fetch predictions'	✓ PASS

Test Cases - Frontend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
1	Renders all form fields	1. Render <LoginForm onLoginSuccess={fn} />. 2. Query for email input, password input, Login button.	Email, password inputs and Login button present	✓ PASS
2	Empty email shows error	1. Click the Login button without filling any field.	Submitting with no email shows 'Email is required'	✓ PASS
3	Empty password shows error	1. Type a valid email. 2. Click Login without filling password.	Submitting with no password shows 'Password is required'	✓ PASS
4	No fetch on validation failure	1. Click Login with empty fields.	fetch() not called when form is invalid	✓ PASS
5	Fetch called with correct payload	1. Type 'User@Example.COM' in email field. 2. Type 'password123' in password field. 3. Click Login.	Email lowercased, sent to /api/v1/auth/login	✓ PASS
6	Token stored in localStorage	1. Submit valid credentials. 2. Check localStorage after submission.	localStorage.getItem('token') equals returned JWT	✓ PASS
7	onLoginSuccess called on success	1. Submit valid credentials. 2. Wait for async resolution.	Callback fires once after successful login	✓ PASS
8	Server error shown	1. Submit credentials. 2. Wait for response.	'Invalid email or password' shown on 401 response	✓ PASS
9	Network error shown	1. Submit valid credentials.	'Network Error' shown when fetch throws	✓ PASS
10	Button disabled while loading	1. Submit valid credentials. 2. Immediately inspect button state.	Submit button is disabled while fetch is pending	✓ PASS
11	Renders all form fields	1. Render <RegisterForm onRegistrationSuccess={fn} />. 2. Query for email, password, confirm password, Register button.	Email, password, confirm password, Register button present	✓ PASS
12	Empty email shows error	1. Click Register without filling any field.	'Email is required' shown on empty submit	✓ PASS
13	Short password shows error	1. Type email. 2. Type password shorter than 8 characters. 3. Click Register.	'at least 8 characters' shown for short password	✓ PASS
14	Password mismatch shows error	1. Fill email and password. 2. Type different value in confirm password. 3. Click Register.	'Passwords do not match' shown when fields differ	✓ PASS
15	No fetch on validation failure	1. Click Register with empty form.	fetch() not called when form is invalid	✓ PASS

Test Cases - Frontend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
16	Fetch called with correct payload	1. Fill 'User@Example.COM', 'password123', 'password123'. 2. Click Register.	Email lowercased, sent to /api/v1/auth/register	✓ PASS
17	Success message and redirect	1. Submit valid form. 2. Advance timers by 800ms.	'Registration successful' shown; onRegistrationSuccess fires after 800ms	✓ PASS
18	Form cleared on success	1. Submit valid form. 2. Wait for success state.	All fields empty after successful registration	✓ PASS
19	Server error shown	1. Submit valid form.	'User already exists' shown on 409 response	✓ PASS
20	Button disabled while submitting	1. Submit valid form. 2. Inspect button state.	Register button disabled while fetch pending	✓ PASS
21	Empty state message shown	1. Render <ImageGallery images={} />.	'No images to display' shown for empty array	✓ PASS
22	Empty state for undefined prop	1. Render <ImageGallery images={undefined} />.	'No images to display' shown when images is undefined	✓ PASS
23	Renders img per image	1. Render <ImageGallery images={[img1, img2]} />. 2. Query all elements.	One element rendered per image with url	✓ PASS
24	title used as alt text	1. Render gallery with image { url: '...', title: 'Knee X-ray' }.	img alt equals image.title	✓ PASS
25	Falls back to originalName	1. Render gallery with image { url: '...', originalName: 'scan.png' }.	img alt equals originalName when title absent	✓ PASS
26	Falls back to fileUrl	1. Render gallery with image { fileUrl: '/uploads/scan.png', title: 'Scan' }.	img src equals fileUrl when url absent	✓ PASS
27	No preview placeholder	1. Render gallery with image { id: '1', title: 'Mystery' }.	'No preview' shown when url and fileUrl both absent	✓ PASS
28	Label text rendered below image	1. Render gallery with image { url: '...', title: 'Knee Scan' }.	Title/name text visible below the image	✓ PASS
29	Renders dropzone text	1. Render <ImageUploader />.	'Drag & drop or click' text visible	✓ PASS
30	Shows drag-active text	1. Fire dragOver event on the dropzone element.	'Drop files here' shown on dragenter	✓ PASS

Test Cases - Frontend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
31	No upload button initially	1. Render <ImageUploader /> without dropping files.	Upload button absent before files are dropped	✓ PASS
32	Upload button after drop	1. Fire drop event with a PNG file. 2. Query for Upload button.	Upload button appears after files are dropped	✓ PASS
33	DICOM placeholder shown	1. Fire drop event with a .dcm file. 2. Query for DICOM label.	'DICOM' placeholder shown for .dcm files	✓ PASS
34	Full upload flow succeeds	1. Drop a PNG file. 2. Click Upload. 3. Wait for all three fetch calls.	presign → PUT → register called in order; onUploadSuccess fires	✓ PASS
35	Alert on presign failure	1. Drop a file. 2. Click Upload.	window.alert called when presign returns non-ok	✓ PASS
36	Button disabled while uploading	1. Drop a file. 2. Click Upload. 3. Inspect button state.	Upload button disabled while fetch is pending	✓ PASS
37	LoginPage — brand heading	1. Render <LoginPage onLoginSuccess={fn} onSwitchToRegister={fn} />.	'KneeVision' h1 rendered	✓ PASS
38	LoginPage — subtitle	1. Render LoginPage.	'Sign in to your account' text present	✓ PASS
39	LoginPage — form rendered	1. Render LoginPage.	Email and password inputs present inside page	✓ PASS
40	LoginPage — switch to register	1. Click the Register button.	onSwitchToRegister called when Register button clicked	✓ PASS
41	RegisterPage — brand heading	1. Render <RegisterPage onRegistrationSuccess={fn} onSwitchToLogin={fn} />.	'KneeVision' h1 rendered	✓ PASS
42	RegisterPage — subtitle	1. Render RegisterPage.	'Create your account' text present	✓ PASS
43	RegisterPage — form rendered	1. Render RegisterPage.	Email and confirm password inputs present inside page	✓ PASS
44	RegisterPage — switch to login	1. Click the Login button.	onSwitchToLogin called when Login button clicked	✓ PASS
45	UploadPage — brand heading	1. Render <UploadPage onLogout={fn} />.	'KneeVision' h1 rendered	✓ PASS


Test Cases - Frontend

Test Case	Test Name	Test Steps	Description / Expected Result	Status
46	UploadPage — subtitle	1. Render UploadPage.	'Upload and review knee X-ray' text present	✓ PASS
47	UploadPage — logout button	1. Render UploadPage.	Logout button visible	✓ PASS
48	UploadPage — logout callback	1. Click the Logout button.	onLogout called when Logout clicked	✓ PASS
49	UploadPage — gallery heading	1. Render UploadPage.	'Gallery' heading present	✓ PASS
50	UploadPage — empty state	1. Render UploadPage with no images.	'No images to display yet' shown initially	✓ PASS
51	UploadPage — add mock image	1. Click 'Add Mock Image'. 2. Query for img element.	Mock image appears in gallery after clicking 'Add Mock Image'	✓ PASS
52	UploadPage — empty state hidden	1. Click 'Add Mock Image'. 2. Query for empty state text.	Empty state message gone after image added	✓ PASS

Bug Fixes

Bugs ID	Bug Name	Description / Expected Result	Status
1	Image Not Found	Uploaded image not getting displayed after session Completion	✓ Fixed
2	Layout discrepancy	Overlapping Widgets on Other Components	✓ Fixed
3	Bread Crumbs in Navigation	Uneven button arrangement for show/hide	✓ Fixed
4	Theme Colour	Text getting blurred while change Dark theme	✓ Fixed
5	Dark Mode in Quiz Page	Did Support Properly	✓ Fixed
6	Navigation Bar	It Need to sticky to position/ Not an Expected Behavior	✓ Fixed



Metrics - Team Velocity

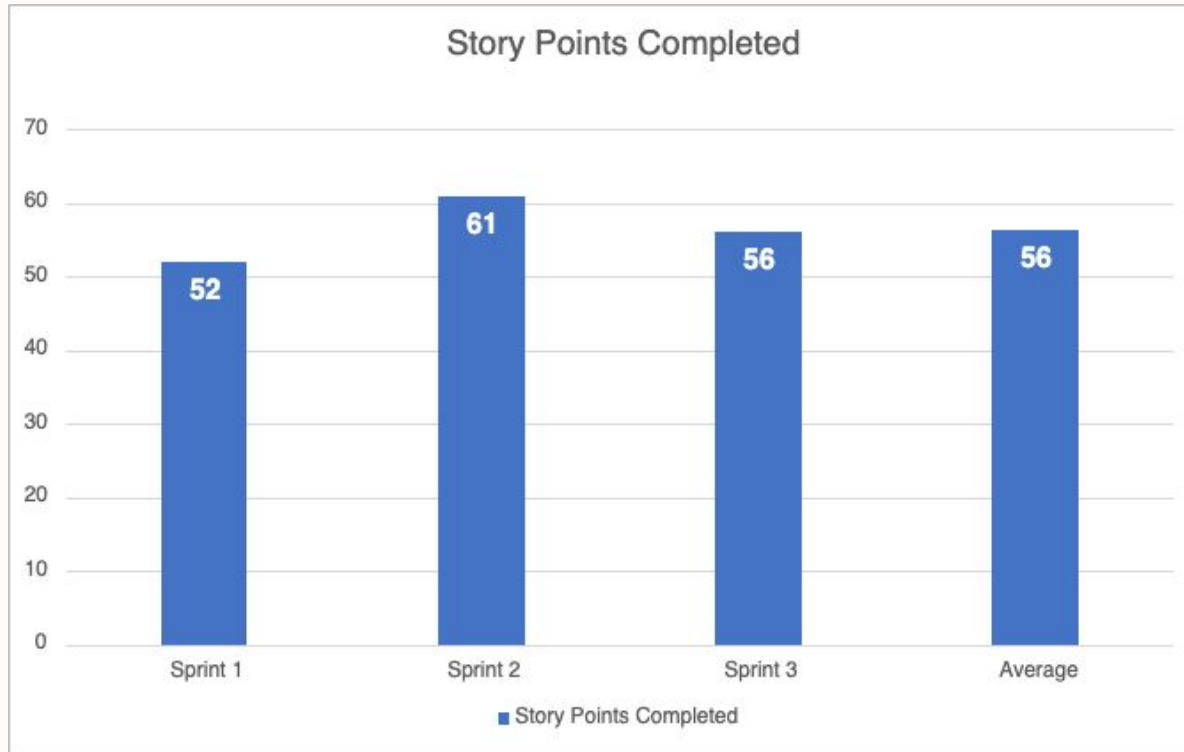



Sprint 3 Story Points Completed: **56**

Average Story Points Completed: **56**


Metrics - Burndown Chart

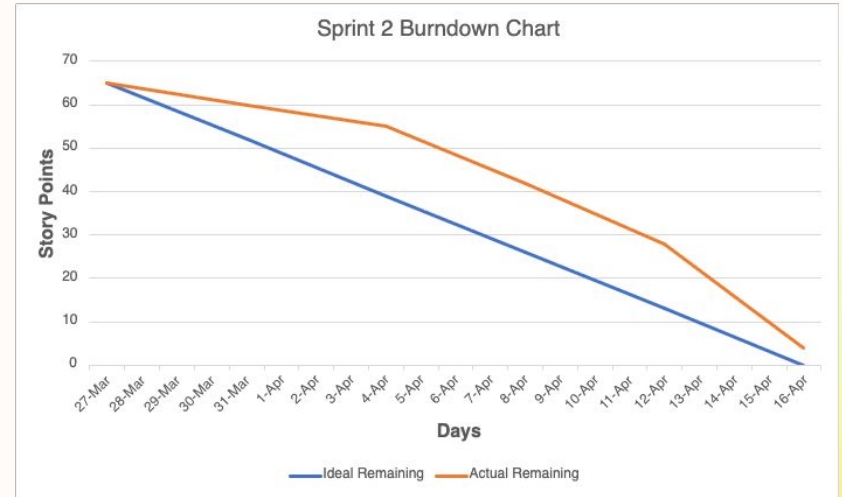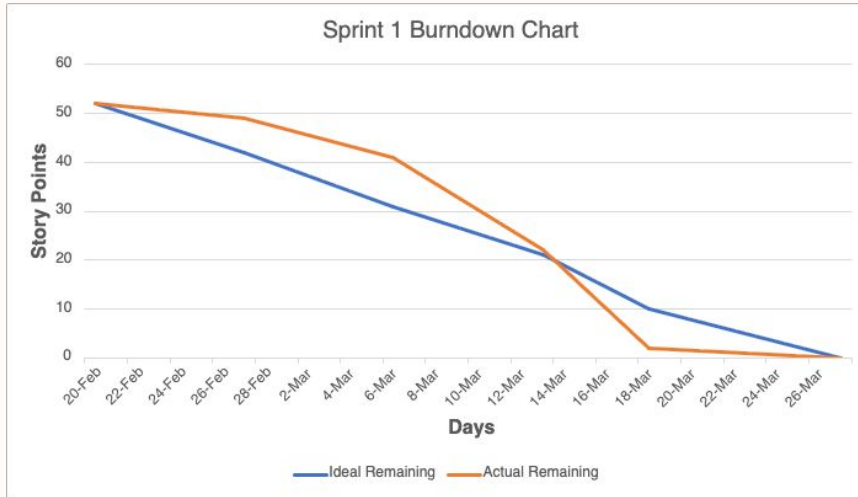

Metrics - Burndown Chart

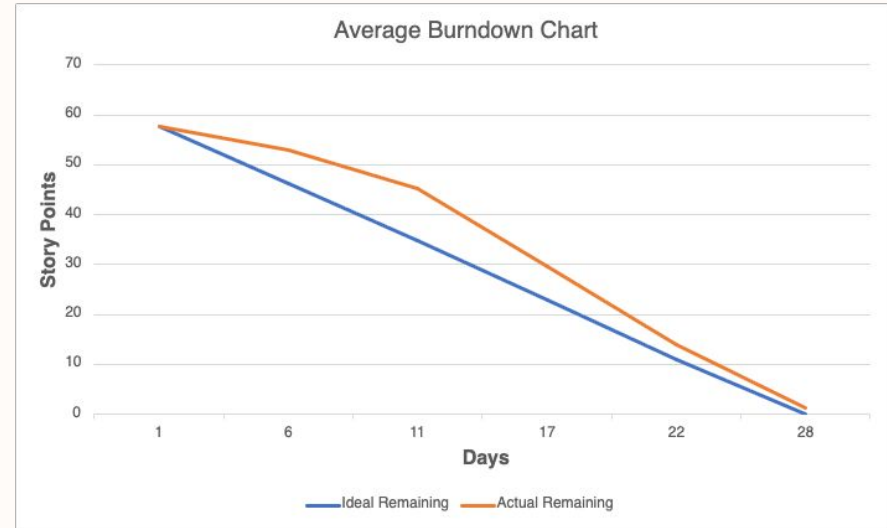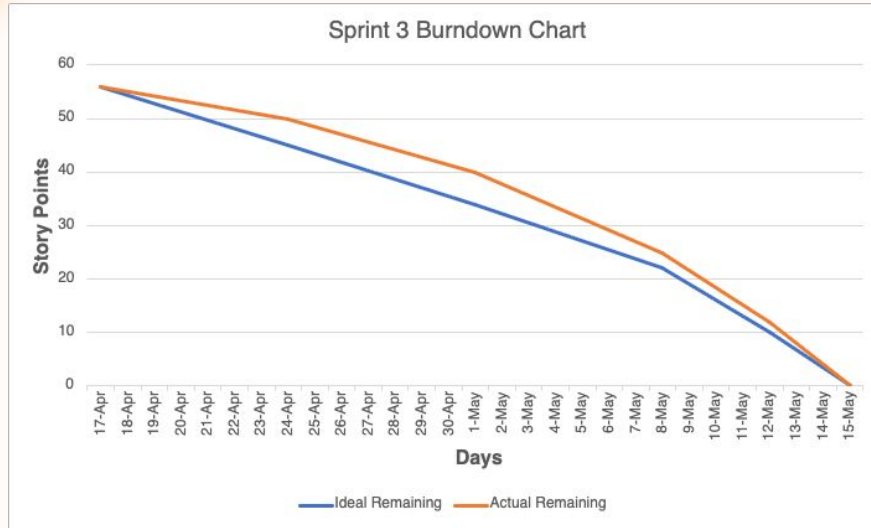

Metrics - Completed/Committed Ratio

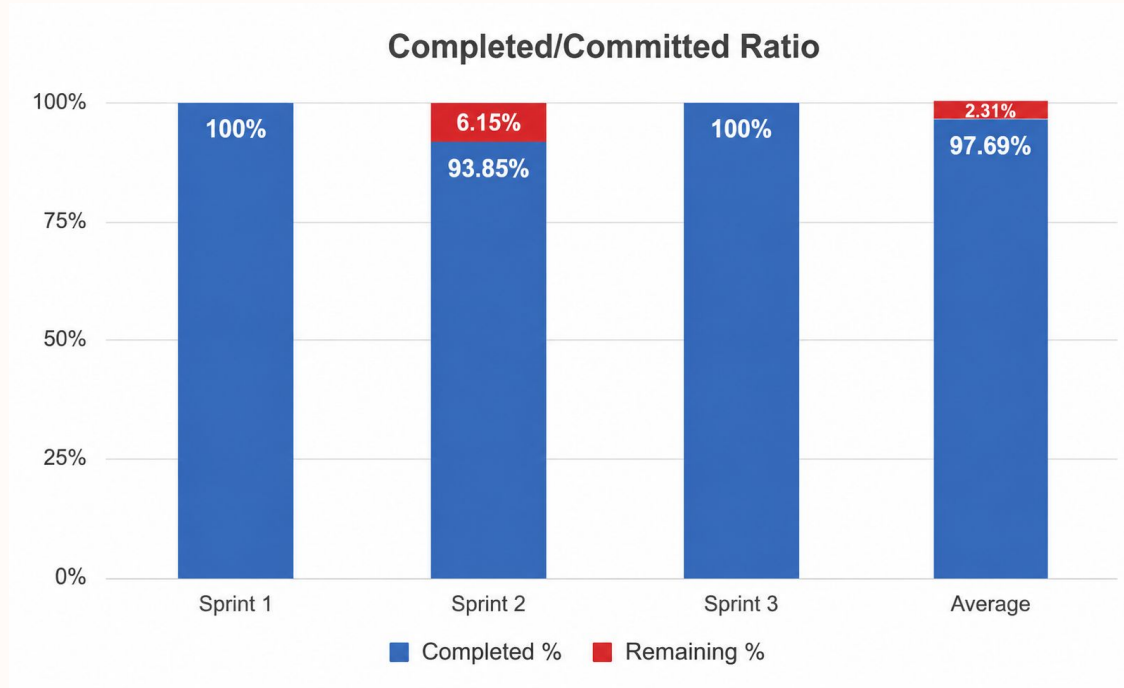

Sprint 3 Story Points Completed: **56/56 (100%)** Average Story Points Completed: **169/173 (97.69%)**

Sprint 3 Retrospective

KneeVision - Sprint 3 Retrospective

What went well +

Base64 image persistence successfully resolved file loss issues in deployment (Render)

+3

End-to-end workflow stabilized (upload → predict → save → retrieve → chat)

+6

UI/UX significantly improved (Profile page, landing page, layout consistency)

+6

Strong integration across services (ML, backend, frontend, chat features)

+4

Team collaboration improved with smoother PR reviews and fewer merge conflicts

+5

Application became demo-ready with all core features functioning reliably

+2

What can be improved +

Prediction pipeline still depends on local /uploads (not fully cloud-based)

+3

Some features (profile updates, password changes) not fully connected to backend

+3

Inconsistent handling of image data (base64 vs file system) across features

+2

Limited time for deeper testing of edge cases and error handling

+5

UI consistency improved but could be further standardized across all pages

+4

Some APIs are tightly coupled and could benefit from cleaner separation of concerns

+4

Action Items +

Refactor prediction pipeline to support full cloud storage (remove /uploads dependency)

+2

Fully connect profile features (username, password, photo) to backend APIs

+2

Standardize data handling strategy for images (base64 vs file-based approach)

+6

Increase test coverage, especially for edge cases and failure scenarios

+4

Continue improving UI consistency using shared styles/components

+1

Refactor APIs for better modularity and maintainability

+4

Sprint 3 Retrospective: Important Points



What Went Well

- End-to-end workflow stabilized (upload → predict → save → retrieve → chat)
- UI/UX significantly improved (Profile page, landing page, layout consistency)
- Team collaboration improved with smoother PR reviews and fewer merge conflicts



Improvements

- Limited time for deeper testing of edge cases and error handling
- UI consistency improved but could be further standardized across all pages
- Some APIs are tightly coupled and could benefit from cleaner separation of concerns



Action Items

- Standardize data handling strategy for images (base64 vs file-based approach)
- Increase test coverage, especially for edge cases and failure scenarios
- Refactor APIs for better modularity and maintainability

Project Roadmap

Sprint Schedule & Key Deliverable

0

DEADLINE: FEB 19

Foundation

Project Charter & MVP
GitHub & Tech Stack Setup
Backlog Creation

DELIVERABLES:

Presentation + Wiki + Backlog

1

DEADLINE: MAR 26

Core Features

User Auth
Image Upload
Image Gallery

DELIVERABLES:

Working Demo + Wiki + Tests

2

DEADLINE: APR 16

AI Engine

Model Training
Prediction API
KL Grade Results

DELIVERABLES:

AI Demo + Draft Paper + Wiki

3

DEADLINE: MAY 14

Deployment

History & CSV Export
Security & UAT Audit
Production Launch

DELIVERABLES:

Final MVP + Full Documentation

Project - APIs - Predict API

```
app.post("/api/v1/predict", requireAuth, async (req, res) => {
  const { fileUrl } = req.body;

  const filename = path.basename(fileUrl);
  const localFilePath = path.join("uploads", filename);

  const prediction = await callMLPredict(localFilePath);

  const imageDoc = await db.collection("images").findOne({ fileUrl });

  const existing = await db.collection("predictions").findOne({
    imageId: imageDoc._id.toString(),
    userId: req.user.userId,
  });

  if (existing) return res.json(existing.result);

  const saved = await db.collection("predictions").insertOne({
    userId: req.user.userId,
    imageId: imageDoc._id.toString(),
    result: prediction,
  });

  res.json({ id: saved.insertedId, ...prediction });
});
```

Project - APIs - Similar Cases API

```
app.post("/api/v1/similar", requireAuth, async (req, res) => {
  const { fileUrl, grade } = req.body;

  const filename = path.basename(fileUrl);
  const fileBuffer = await fs.readFile(path.join("uploads", filename));

  const form = new FormData();
  form.append("file", new Blob([fileBuffer]), filename);

  const response = await fetch(`${ML_SERVICE_URL}/similar`, {
    method: "POST",
    body: form,
  });

  const data = await response.json();
  res.json(data);
});
```


Project - APIs - Chat API

```
app.post("/api/v1/chat", requireAuth, async (req, res) => {
  const { imageBase64, result, messages } = req.body;

  const systemPrompt = `KL Grade: ${result.grade}
Confidence: ${result.confidence}`;

  const chatMessages = [
    { role: "system", content: systemPrompt },
    {
      role: "user",
      content: [
        { type: "image_url", image_url: { url: imageBase64 } },
        { type: "text", text: messages[0].content },
      ],
    },
  ];

  const response = await fetch("https://api.groq.com/openai/v1/chat/completions", {
    method: "POST",
    headers: { Authorization: `Bearer ${process.env.GROQ_API_KEY}` },
    body: JSON.stringify({
      model: "llama-4",
      messages: chatMessages,
    }),
  });

  const data = await response.json();
  res.json({ reply: data.choices[0].message.content });
});
```

Docker:

Containers [Give feedback](#)

View all your running containers and applications. [Learn more](#)

Container CPU usage ⓘ

0.64% / 1200% (12 CPUs available)

Container memory usage ⓘ

606.64MB / 7.44GB

[Show charts](#)


Q Search




Only show running containers

☐	Name	Container ID	Image	Port(s)	CPU (%)	Last started	Actions
☐ ▾ ●	infra	-	-	-	0.64%	43 minutes ago	
☐ ●	ml-service-1	90d81254ea8b	infra-ml-service	8000:8000 ↗	0.34%	43 minutes ago	
☐ ●	backend-1	81a2925d5724	infra-backend	4000:4000 ↗	0.3%	43 minutes ago	
☐ ●	frontend-1	d0a7fbdeffc1	infra-frontend	80:80 ↗	0%	43 minutes ago	

Project Demo - Application Screenshots

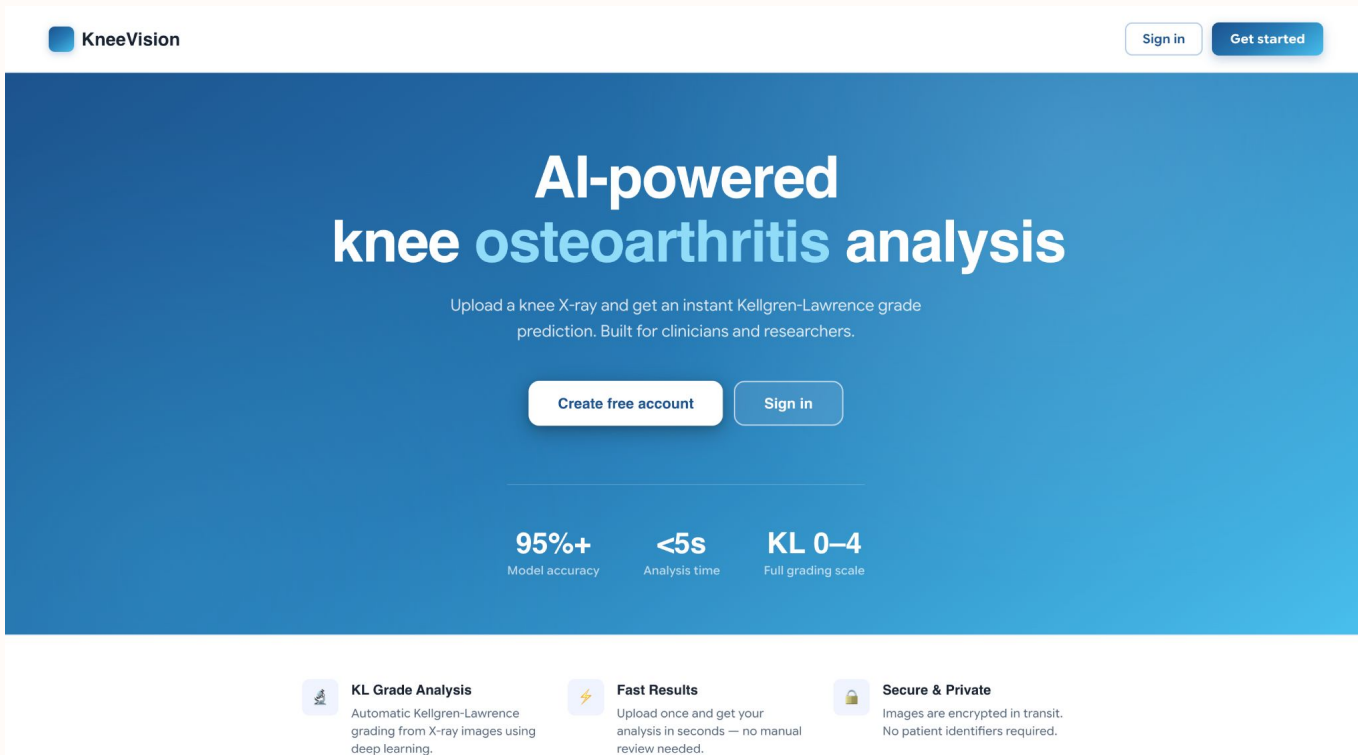

Project Demo - Application Screenshots

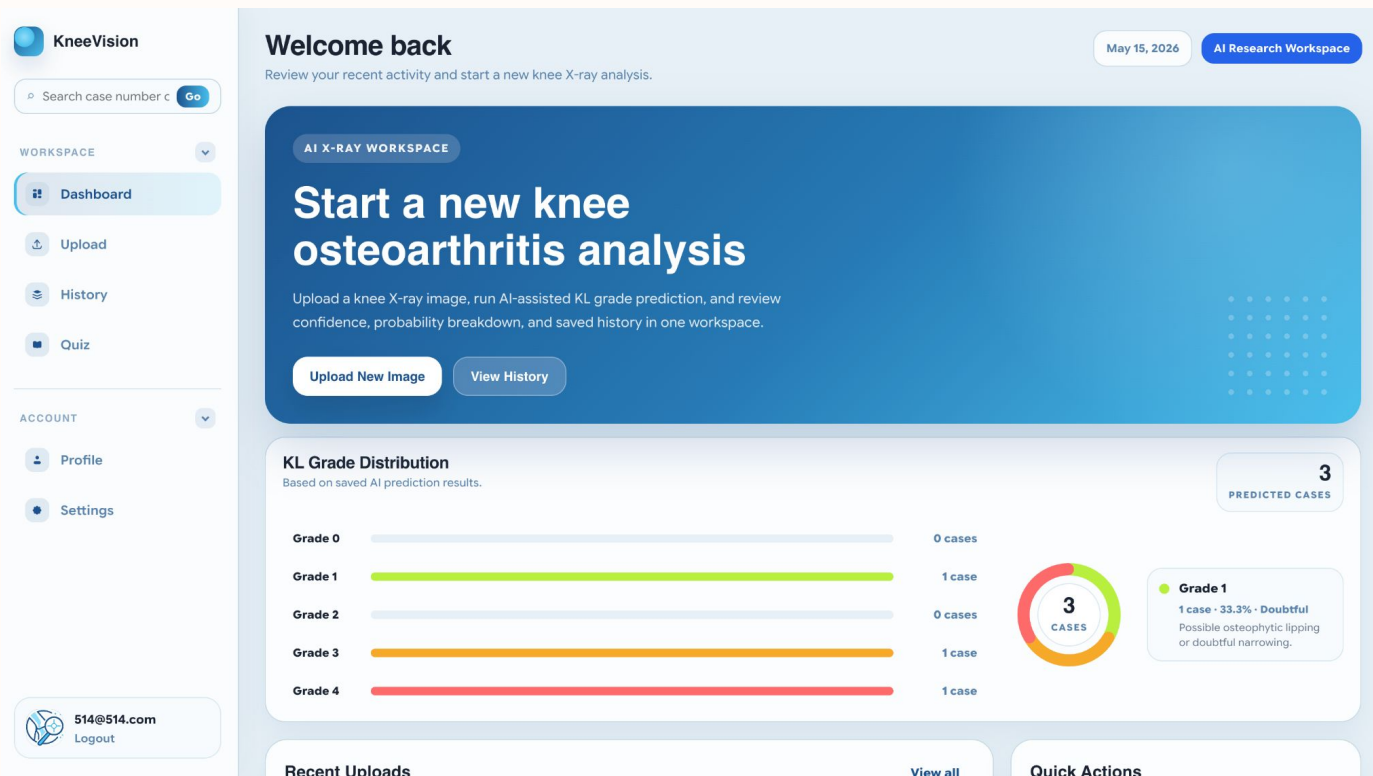

Project Demo - Application Screenshots

 **KneeVision**

Go

WORKSPACE 

 Dashboard

 Upload

 History

 Quiz

ACCOUNT 

 Profile

 Settings

 **514@514.com**
Logout

Upload X-ray

Upload a knee X-ray image and run AI-assisted osteoarthritis severity analysis.

New X-ray Analysis

Select a JPG or PNG knee X-ray image. After upload, KneeVision will automatically run the AI model and open the result page.




Drag & drop or click to browse
Supports JPG, PNG, WebP, DICOM · batch uploads OK

JPG

PNG

WebP

DICOM

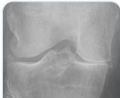

9814262R.png

Upload

1 file selected

Gallery

Browse, filter, and reopen previous uploads.

FILTER 

All 

SORT 

Newest first 

May 15, 2026

AI Research Workspace

Upload Guide

Recommended input checklist for better analysis quality.

01

Use clear X-ray images
Avoid heavily blurred or cropped knee joint images.


02

Review confidence score
Confidence supports interpretation, but it is not diagnosis.

03

Check saved history
Uploaded cases and predictions can be reviewed later.

Project Demo - Application Screenshots

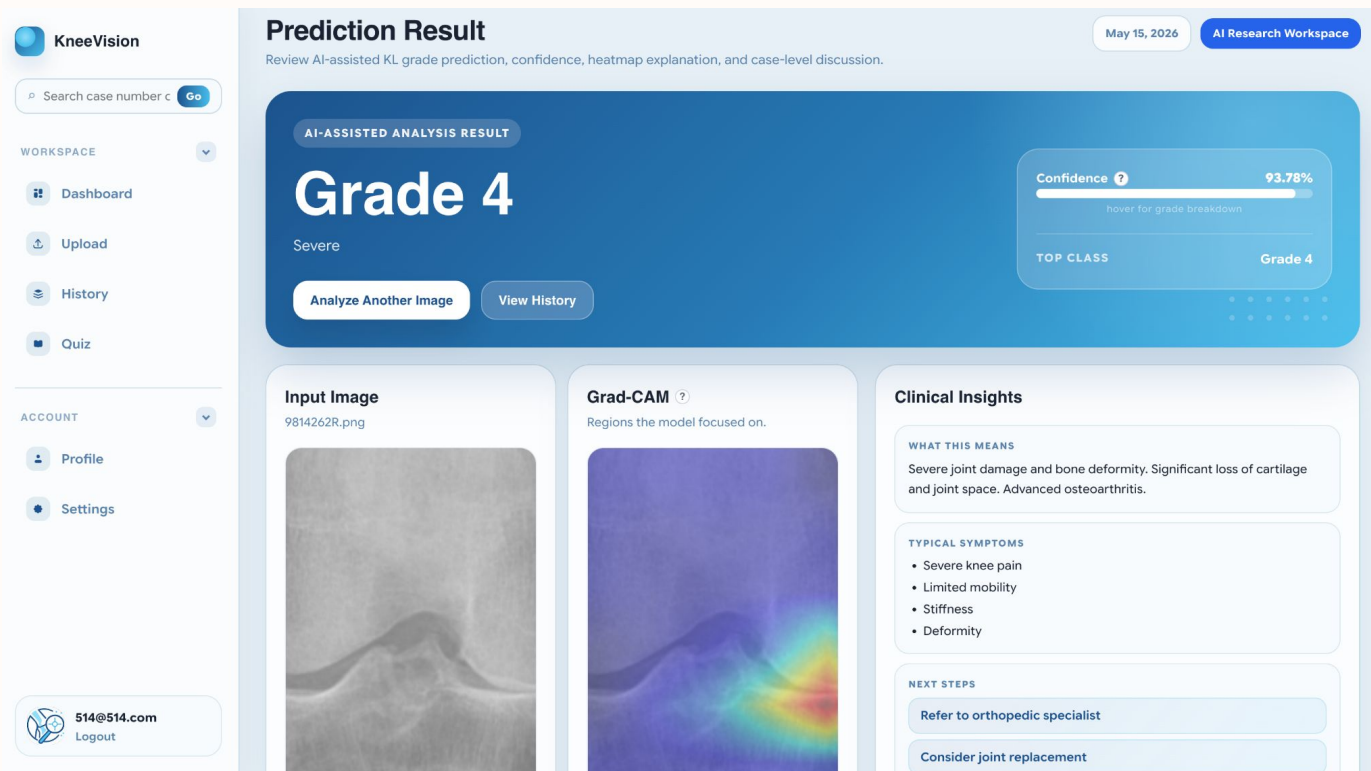

Project Demo - Application Screenshots

 **KneeVision**

Go

WORKSPACE

 Dashboard

 Upload

 History

 Quiz

ACCOUNT


 Profile

 Settings

 **514@514.com**
Logout

Similar Cases 

Top 5 visually similar X-rays from the reference database. Click a case for clinical details.

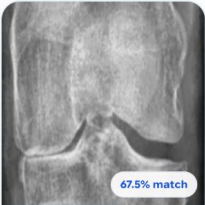


67.5% match

KL Grade 4 · Grade 4

Risk: High risk of progression to end-stage disease

65-year-old male

Click to view details

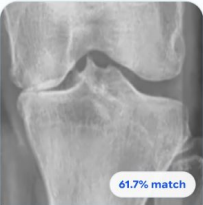


61.7% match

KL Grade 4 · Grade 4

Risk: High risk of further joint damage

Female

Click to view details

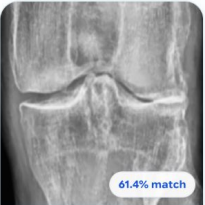


61.4% match

KL Grade 4 · Grade 4

Risk: High risk of further joint damage

Female

Click to view details

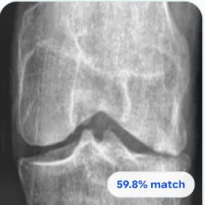


59.8% match

KL Grade 4 · Grade 4

Risk: High risk of progression to joint replacement

Female

Click to view details

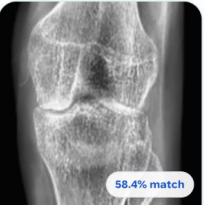

58.4% match

KL Grade 4 · Grade 4

Risk: High

65-year-old female

Click to view details

Ask about this X-ray 

Ask questions about the predicted grade, findings, or similar cases.

Try: "Why is this Grade 2?" or "How does this compare to the similar cases?"

Project Demo - Application Screenshots

The screenshot displays the KneeVision application interface. On the left is a sidebar with navigation options: Dashboard, Upload, History, Quiz, Profile, and Settings. The main content area is divided into two sections. The top section shows a row of five case cards, each titled 'KL Grade 4 · Grade 4' and containing a risk assessment, gender, and age. The bottom section is titled 'Ask about this X-ray' and contains a detailed AI analysis of an X-ray image, including a description of the findings, a list of characteristics for Grade 4, and a confidence score.

KneeVision

Search case number c **Go**

WORKSPACE

- Dashboard
- Upload
- History
- Quiz

ACCOUNT

- Profile
- Settings

KL Grade 4 · Grade 4
Risk: High risk of further joint damage
Female
Click to view details

KL Grade 4 · Grade 4
Risk: High risk of progression to end stage
65-year-old male
Click to view details

KL Grade 4 · Grade 4
Risk: High risk of further joint damage
Female
Click to view details

KL Grade 4 · Grade 4
Risk: High risk of progression to joint replacement
Female
Click to view details

KL Grade 4 · Grade 4
Risk: High risk of further joint damage
Female
Click to view details

Ask about this X-ray ?
Ask questions about the predicted grade, findings, or similar cases. · [History restored](#)

Why is this Grade 4?

KNEEVISION AI

The X-ray image shows severe joint space narrowing, which is a hallmark of advanced osteoarthritis. The predicted Grade 4 classification is based on the Kellgren-Lawrence grading system, which evaluates the severity of osteoarthritis based on radiographic features.

Specifically, Grade 4 is characterized by:

- Severe joint space narrowing
- Large osteophytes
- Subchondral sclerosis
- Cystic changes

In this image, I observe severe joint space narrowing, which suggests significant cartilage loss. Although osteophytes, subchondral sclerosis, and cystic changes are not clearly visible in this single X-ray image, the presence of severe joint space narrowing is consistent with a Grade 4 classification.

The model has predicted a Grade 4 classification with 94.61% confidence, indicating a high level of confidence in this assessment.

Ask about this X-ray... **Send**

514@514.com
Logout

Project Demo - Application Screenshots

 **KneeVision**

Go

WORKSPACE 

 Dashboard

 Upload

 History

 Quiz

ACCOUNT 

 Profile

 Settings

 **514@514.com**
Logout

Analysis History

Browse saved uploads, filter AI predictions, and export results for review.

Total Records
4

Predicted Cases
4

Avg Confidence
93.09%

Saved Analysis Records 

Export CSV 

FILTER  All 

GRADE  All Grades 

SORT  Newest first 




9778632R.png
5/15/2026, 8:37:26 AM

Grade 4  90.91% 

Severe 

Delete

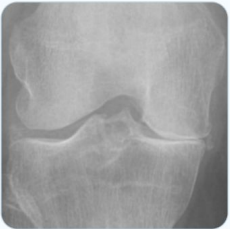


9814262R.png
5/15/2026, 8:13:10 AM

Grade 4  93.78% 

Severe 

Delete

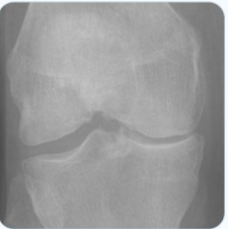

9992318R.png
5/14/2026, 9:07:04 PM

Grade 3  93.04% 

Moderate 

Delete



9924274L.png
5/14/2026, 9:02:29 PM

Grade 4  94.61% 

Severe 

Delete

Project Demo - Application Screenshots

The screenshot displays the 'KneeVision Quiz' application interface. On the left is a sidebar with the 'KneeVision' logo, a search bar for case numbers, and navigation links for 'Workspace' (Dashboard, Upload, History, Quiz) and 'Account' (Profile, Settings). The main content area is titled 'KneeVision Quiz' with a subtitle 'Practice Kellgren-Lawrence grading using real knee X-ray cases.' It includes a date 'May 15, 2026' and a 'AI Research Workspace' badge. A large blue banner for the 'KL GRADING CHALLENGE' features the text 'Train your eye for knee OA severity' and instructions to review X-ray images and choose the most likely KL grade. Below this are 'Start Quiz' and 'View Leaderboard' buttons. The 'Choose Difficulty' section offers three options: 'Easy' (Grade 0 vs Grade 4 only), 'Medium' (All 5 grades), and 'Hard' (Grades 1, 2, and 3 only). The 'Quiz Rules' section lists four rules: 10 questions per session, 30 seconds per image, building streaks for consecutive correct answers, and saving top scores to the leaderboard.

KneeVision

Search case number c **Go**

WORKSPACE

- Dashboard
- Upload
- History
- Quiz**

ACCOUNT

- Profile
- Settings

KneeVision Quiz

Practice Kellgren-Lawrence grading using real knee X-ray cases.

May 15, 2026 **AI Research Workspace**

KL GRADING CHALLENGE

Train your eye for knee OA severity

Review knee X-ray images, choose the most likely KL grade, and get immediate feedback based on radiographic findings.

Start Quiz **View Leaderboard**

Choose Difficulty

Select the grade range you want to practice.

- Easy**
Grade 0 vs Grade 4 only
- Medium**
All 5 grades
- Hard**
Grades 1, 2, and 3 only

Quiz Rules

Each session is short and focused.

- 01** 10 questions per session
- 02** 30 seconds per image
- 03** Build streaks for consecutive correct answers
- 04** Top scores are saved to the leaderboard

514@514.com
Logout

Project Demo - Application Screenshots

 **KneeVision**

Go

WORKSPACE 

 Dashboard

 Upload

 History

 Quiz

ACCOUNT 

 Profile

 Settings

 **514@514.com**
Logout

KneeVision Quiz

Question 1 of 10. Choose the most likely KL grade before time runs out.

Question 1 / 10

Score: 0

3s remaining



What is the KL Grade?

Select the best answer for this X-ray image.

Grade 0 — Normal

Grade 1 — Doubtful

Grade 2 — Mild

Grade 3 — Moderate

Grade 4 — Severe

Project Demo - Application Screenshots

The screenshot displays the 'KneeVision Quiz' interface. On the left is a sidebar with a 'KneeVision' logo, a search bar, and navigation links for 'Workspace' (Dashboard, Upload, History, Quiz) and 'Account' (Profile, Settings). The main area shows a quiz result for 'Grade 0 — Normal' with a 'Correct!' message. Below this, patient information is presented in a structured layout, including progression risk, recommended follow-up, patient profile, lifestyle factors, previous findings, and actions. A score of 1/10 and a 'Next Question' button are also visible.

KneeVision

Search case number c **Go**

WORKSPACE

- Dashboard
- Upload
- History
- Quiz**

ACCOUNT

- Profile
- Settings

KneeVision Quiz

Review the explanation, then move to the next case.

May 15, 2026 **AI Research Workspace**

Correct!

Correct Answer: Grade 0 — Normal

Review the radiographic findings for this case.

PROGRESSION RISK Low	RECOMMENDED FOLLOW-UP Annual X-ray in 5 years
PATIENT PROFILE 45-year-old female, 65 kg, 170 cm, no significant medical history	LIFESTYLE FACTORS Good diet, regular exercise, normal weight
Previous: None	
Actions: Maintain current lifestyle and exercise routine	

Score: 1 Question 3 of 10

Next Question →

514@514.com
Logout

Project Demo - Application Screenshots

KneeVision

Go

WORKSPACE

Dashboard

Upload

History

Quiz

ACCOUNT

Profile

Settings

514@514.com
Logout

Profile

Manage your KneeVision account details and security.

Account Information

Your email address and display name on this device.

Email

514@514.com

Username

Not set yet — choose a display name

Change Username

Save Username

Profile Photo

Shown in the sidebar. Stored locally in your browser.



Upload Photo

Remove Photo

JPG or PNG recommended. Stored locally for this MVP.

Security

Update your account password.

Current Password

New Password

Confirm New Password

Session & Security

Your current login session details.

Logged In As

514@514.com

Session Status

Active

Last Login

May 15, 2026, 8:06 AM

May 15, 2026

AI Research Workspace

Project Demo - Application Screenshots

The screenshot displays the 'Settings' page of the KneeVision application. The interface is divided into a left sidebar and a main content area. The sidebar includes a search bar, a 'Workspace' section with links to Dashboard, Upload, History, and Quiz, and an 'Account' section with links to Profile and Settings. The main content area is titled 'Settings' and includes a subtitle 'Adjust local display preferences for your KneeVision workspace.' It features four primary settings panels: 'Appearance' (with 'Light Mode' selected), 'Interface Preferences' (with three toggle switches), 'Language' (showing 'Current Language' as English), and 'System Info' (displaying version, frontend, backend, and ML service details). A 'Reset Local Preferences' button is located at the bottom of the settings area. The top right of the page shows the date 'May 15, 2026' and the workspace name 'AI Research Workspace'.

KneeVision

Search case number c

WORKSPACE

- Dashboard
- Upload
- History
- Quiz

ACCOUNT

- Profile
- Settings

514@514.com Logout

Settings

Adjust local display preferences for your KneeVision workspace.

Appearance

Choose how the workspace should look on this device.

Light Mode
Clean blue medical workspace

Dark Mode
Low-light review interface

Interface Preferences

These preferences are saved locally in this browser.

- Educational Tips**
Show short explanations for medical and AI terms. ☒
- Compact Layout**
Reduce spacing slightly for smaller screens. ☐
- Remember Last Page**
Keep this as a local preference for future navigation. ☐

Language

Language support can be expanded in a future version.

Current Language
English

System Info

Basic deployment information for this project build.

Version v1.0.0	Frontend Vercel
Backend Render	ML Service Hugging Face

Reset Local Preferences

Project Demo - Application Screenshots

The screenshot displays the 'KneeVision' AI Research Workspace. The interface is divided into a left sidebar, a top header, and a main content area.

Left Sidebar:

- KneeVision** logo and name.
- Search bar: "Search case number c" with a "Go" button.
- WORKSPACE** section with a dropdown arrow:
 - Dashboard (with a grid icon)
 - Upload (with an upload icon)
 - History (with a list icon)
 - Quiz (with a document icon)
- ACCOUNT** section with a dropdown arrow:
 - Profile (with a person icon)
 - Settings (with a gear icon)
- Footer: "514@514.com" and "Logout" button.

Top Header:

- AI Research Workspace** title.
- Subtitle: "Explore osteoarthritis research with cited sources from trusted medical databases."
- Date: "May 15, 2026".
- Button: "AI Research Workspace".

Main Content Area:

- Sessions** tab with a "+ New" button. A session titled "Explain the Kellgren-Lawren..." is selected.
- YOU** (User) section with a blue button: "Explain the Kellgren-Lawrence grading system for knee osteoarthritis in detail."
- RESEARCH ASSISTANT** section:
 - SUMMARY**: "The Kellgren-Lawrence (KL) grading system is a widely used radiographic classification system for knee osteoarthritis (OA), assessing the severity of OA based on X-ray features. The system grades knee OA from 0 to 4, with Grade 0 indicating no radiographic features of OA and Grade 4 indicating severe OA. This grading system helps clinicians and researchers evaluate the progression of OA and determine the most suitable treatment options."
 - KEY POINTS**:
 - The KL grading system evaluates knee OA based on five key radiographic features: joint space narrowing, osteophytes, subchondral sclerosis, subchondral cysts, and varus/valgus deformity.
 - The system consists of five grades:
 - Grade 0: No radiographic features of OA, normal joint space, and no osteophytes.
 - Grade 1: Doubtful narrowing of joint space, possible osteophytic lipping, but questionable significance.
 - Grade 2: Definite osteophytes, possible narrowing of joint space, indicating mild OA.
 - Grade 3: Multiple osteophytes, definite narrowing of joint space, some sclerosis, possible deformity, indicating moderate OA.
 - Grade 4: Large osteophytes, marked narrowing of joint space, severe sclerosis, definite deformity, indicating severe OA.
 - KL Grade Reference** tab.
- Chat Input**: "Ask about OA grading, treatments, imaging features... (Enter to send, Shift+Enter for new line)" with a "Send" button.

Live Application Demo



Future Enhancements

- **Federated Learning** — train across multiple hospital datasets without sharing patient data, improving model generalization while maintaining privacy compliance (HIPAA/GDPR)
- **Treatment Outcome Prediction** — given a patient's current grade and chosen treatment path, predict the likely outcome at 6/12/24 months using historical case data
- **Mobile Deployment** — quantize the EfficientNet-B4 model (INT8/FP16) for edge deployment on mobile devices, enabling point-of-care screening in rural clinics without internet connectivity

AI Disclosure

- Assisted in generating initial user stories and backlog items
- Helped draft test cases and testing scenarios
- Supported development of specific backend API routes (authentication, image upload, and database integration) and frontend components (forms, state handling, and API calls)
- Provided guidance for debugging and resolving errors
- Assisted with structuring documentation and project artifacts

Wiki Page Link

<https://github.com/htmhw/STARLABS/wiki#kneevision>



THANK YOU

Team S.T.A.R Labs